

OptiFlow Turbo V Ion Source

for the ZenoTOF 7600/7600+ System

Operator Guide



This document is provided to customers who have purchased SCIEX equipment to use in the operation of such SCIEX equipment. This document is copyright protected and any reproduction of this document or any part of this document is strictly prohibited, except as SCIEX may authorize in writing.

Software that may be described in this document is furnished under a license agreement. It is against the law to copy, modify, or distribute the software on any medium, except as specifically allowed in the license agreement. Furthermore, the license agreement may prohibit the software from being disassembled, reverse engineered, or decompiled for any purpose. Warranties are as stated therein.

Portions of this document may make reference to other manufacturers and/or their products, which may contain parts whose names are registered as trademarks and/or function as trademarks of their respective owners. Any such use is intended only to designate such products as those manufacturers' products and does not imply any right and/or license to use or permit others to use such manufacturers' and/or their product names as trademarks.

SCIEX warranties are limited to those express warranties provided at the time of sale or license of its products and are the sole and exclusive representations, warranties, and obligations of SCIEX. SCIEX makes no other warranty of any kind whatsoever, expressed or implied, including without limitation, warranties of merchantability or fitness for a particular purpose, whether arising from a statute or otherwise in law or from a course of dealing or usage of trade, all of which are expressly disclaimed, and assumes no responsibility or contingent liability, including indirect or consequential damages, for any use by the purchaser or for any adverse circumstances arising therefrom.

(GEN-IDV-09-10816-E)

For Research Use Only. Not for use in Diagnostic Procedures.

Trademarks and/or registered trademarks mentioned herein, including associated logos, are the property of AB Sciex Pte. Ltd., or their respective owners, in the United States and/or certain other countries (see sciex.com/trademarks).

AB Sciex™ is being used under license.

Echo, Echo MS, and Echo MS+ are trademarks or registered trademarks of Labcyte, Inc. in the United States and other countries, and are being used under license.

© 2024 DH Tech. Dev. Pte. Ltd.



AB Sciex Pte. Ltd.
Blk33, #04-06 Marsiling Industrial Estate Road 3
Woodlands Central Industrial Estate, Singapore 739256

Contents

1 Operational Precautions and Limitations	5
Operational Precautions and Hazards.....	5
Chemical Precautions.....	6
Laboratory Conditions.....	7
Safe Environmental Conditions.....	7
Performance Specifications.....	7
Equipment Use and Modification.....	8
2 Ion Source Overview	9
Ion Source Components.....	9
Micro Probe.....	9
Nano Probe.....	12
Probe for the Echo [®] MS System.....	13
Gas and Electrical Connections.....	14
Ion Source Sense Circuit.....	14
Source Exhaust System.....	14
3 Ion Source Installation	16
Install the Probe for the Echo [®] MS System or the Micro Probe On the Ion Source.....	17
Install the Ion Source on the Mass Spectrometer.....	19
Install the Electrode in a Micro Probe.....	20
Install an Infusion Adapter and Connect an Infusion Line to a Micro Probe.....	23
Install the Micro Column and Heater.....	25
Install the Electrode in the Nano Probe.....	32
Install the Nano Column Cartridge and Heater.....	35
Connect an Infusion Line to the Nano Probe.....	41
Install ESI Calibration Probe and Electrode.....	43
Install the Calibration Tubing for the Echo [®] MS+ System.....	47
Sample Inlet Requirements.....	50
Examine for Leaks.....	50
4 Ion Source Maintenance	52
Recommended Maintenance Schedule.....	53
OptiFlow Turbo V Ion Source Handling.....	54
Remove the Ion Source.....	57
Clean the Ion Source Surfaces.....	58
Remove the Probe.....	58

Contents

Clean the Electrode	59
Storage and Handling	60
A Connect the Echo[®] MS+ System to the OptiFlow Turbo V Ion Source	61
B Labels on the Ion Source	73
C Glossary of Symbols	74
Contact Us	80
Customer Training	80
Online Learning Center	80
SCIEX Support	80
Cybersecurity	80
Documentation	80

Operational Precautions and Limitations

1

Note: Before operating the system, carefully read all of the sections of this guide.

This section contains general safety-related information. It also describes potential hazards and associated warnings for the system and the precautions that should be taken to minimize the hazards.

For information about the symbols and conventions used in the laboratory environment, on the system, and in this documentation, refer to the section: [Glossary of Symbols](#).

Operational Precautions and Hazards

For regulatory and safety information for the mass spectrometer, refer to the document: *System User Guide*.



WARNING! Ionizing Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Do not use the ion source without knowledge of and training in the proper use, containment, and evacuation of toxic or injurious materials used with the ion source.



WARNING! Hot Surface Hazard. Before any maintenance procedures are started, let the temperature of the OptiFlow Turbo V ion source decrease for at least 60 minutes. Some surfaces of the ion source and vacuum interface become hot during operation.



WARNING! Fire and Toxic Chemical Hazard. Before solvent leaks are corrected, make sure that the liquid flow to the ion source is stopped, that the ion spray voltage is off, that there are no open flames or other sources of fire nearby, and that the room is sufficiently ventilated. The fluid from a leak can be highly flammable. If the fluid is exposed to electrical discharges or a source of fire, then ignition can occur. If ventilation is not sufficient, then the fluid might cause poisoning.



WARNING! Toxic Chemical Hazard. Wear personal protective equipment (PPE), including a laboratory coat, gloves, and safety glasses, to avoid skin or eye exposure.

Operational Precautions and Limitations



WARNING! Ionizing Radiation Hazard, Biohazard, or Toxic Chemical Hazard. In the event of a chemical spill, review product safety data sheets for specific instructions. Make sure that the system is in Standby state before cleaning a spill near the ion source. Use applicable personal protective equipment and absorbent wipes to contain the spill and dispose of it following local regulations.



WARNING! Environmental Hazard. Do not discard system components in municipal waste. To discard components correctly, obey local regulations.



WARNING! Electrical Shock Hazard. Avoid contact with the high voltages applied to the ion source during operation. Put the system in Standby state before adjusting the sample tubing or other equipment near the ion source.

Note: Use zero air when using the OptiFlow Turbo V ion source with micro flow rates under 10 $\mu\text{L}/\text{min}$ or nano flow rates. Do not use UHP nitrogen for Ion source gas 1 or Ion source gas 2, because there is an increased risk of corona discharge, which can damage the emitter tip.

Chemical Precautions



WARNING! Ionizing Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Before cleaning or maintenance, identify whether decontamination is required. If radioactive materials, biological agents, or toxic chemicals have been used with the system, then the customer must decontaminate the system before cleaning or maintenance.



WARNING! Puncture Hazard, Ionizing Radiation Hazard, Biohazard, or Toxic Chemical Hazard. If the ion source window is cracked or broken, then do not use the ion source. Contact a SCIEX field service employee (FSE). Any toxic or injurious materials introduced in the equipment will be present in the source exhaust output. Exhaust from equipment should be vented from the room. Dispose of sharps following established laboratory safety procedures.



WARNING! Environmental Hazard. Do not discard system components in municipal waste. To discard components correctly, obey local regulations.



WARNING! Biohazard or Toxic Chemical Hazard. To prevent leaks, connect the drain tubing to the mass spectrometer and the source exhaust drain bottle correctly.

- Before servicing and regular maintenance, identify the chemicals that have been used in the system. For the health and safety precautions that must be obeyed for a chemical, refer to the

safety data sheet (SDS). For storage information, refer to the certificate of analysis. To find a SCIEX SDS or certificate of analysis, go to [sciex.com/tech-regulatory](https://www.sciex.com/tech-regulatory).

- Always wear assigned personal protective equipment, which includes powder-free gloves, protective eyewear, and a laboratory coat.

Note: Nitrile or neoprene gloves are recommended.

- Do work in a well-ventilated area or fume hood.
- When flammable materials such as isopropanol and methanol are in use, do not go near ignition sources.
- Be careful with the use and disposal of any chemicals. If the correct procedures for chemical use and disposal are not obeyed, then personal injury can occur.
- During cleaning, do not let chemicals touch the skin. Wash hands after use.
- Make sure that all exhaust hoses are connected correctly and that all connections are functioning as designed.
- Collect all spent liquids and discard them as hazardous waste.
- Obey all of the local regulations for the storage, use, and disposal of biohazardous, toxic, and radioactive materials.

Laboratory Conditions

Safe Environmental Conditions

The system is designed to operate safely under these conditions:

- Indoors
- Altitude: Up to 2,000 m (6,560 ft) above sea level
- Ambient temperature: 10 °C (50 °F) to 35 °C (95 °F)
- Relative humidity: 20% to 80%, non-condensing
- Mains supply voltage fluctuations: $\pm 10\%$ of the nominal voltage
- Transient overvoltages: Up to the levels of Overvoltage Category II
- Temporary overvoltages on the mains supply
- Pollution Degree 2

Performance Specifications

The system is designed to meet specifications under these conditions:

- Ambient temperature is between 15 °C to 30 °C (59 °F to 86 °F).

Operational Precautions and Limitations

Over time, the temperature must remain within a range of 2 °C (3.6 °F), with the rate of the change in temperature not exceeding 2 °C (3.6 °F) per hour. Ambient temperature fluctuations exceeding the limits might result in mass shifts in spectra.

- Relative humidity is between 20% to 80%, non-condensing.

Equipment Use and Modification



WARNING! Electrical Shock Hazard. Do not remove the covers. If the covers are removed, then injury or incorrect system operation can occur. Removal of the covers is not required for routine maintenance, inspection, or adjustment. For repairs that require removal of the covers, contact a SCIEX field service employee (FSE).



WARNING! Personal Injury Hazard. Use only parts that are recommended by SCIEX. The use of parts that are not recommended by SCIEX or the use of parts for any purpose other than their intended purpose can put the user at risk of harm or have a negative effect on system performance.

Use the system indoors in a laboratory that complies with the environmental conditions recommended in the mass spectrometer document: *Site Planning Guide*.

If the system is used in conditions or in an environment that are not approved by the manufacturer, then the performance and protection that is supplied by the equipment might be decreased or lost.

Contact an FSE for information about servicing the system. Unauthorized modification or operation of the system might cause personal injury and equipment damage, and might void the warranty. If the system is operated outside the recommended environmental conditions or with unauthorized modifications, then the acquired data might be inaccurate.

Ion Source Overview

2

The OptiFlow Turbo V ion source can be used with the ZenoTOF 7600 or ZenoTOF 7600+ system.

Note: The ZenoTOF 7600+ system does not support the Echo[®] MS+ system.

The ion source has two probe ports, a front port and a top port.

- If the ion source is configured for micro flow operation, then a Micro probe must be installed in the top port, and the ESI calibration probe or a probe port plug must be installed in the front port. Refer to the figure: [Figure 2-1](#).
- If the ion source is configured for nano flow operation, then a Nano probe must be installed in the front port, and a probe port plug must be installed in the top port. Refer to the figure: [Figure 2-2](#).

The probes are configured for various flow rates and electrodes:

- Micro 1–50 μL probe: Flow rate of 1 $\mu\text{L}/\text{min}$ to 50 $\mu\text{L}/\text{min}$. Compatible electrodes are:
 - Electrode 1–10 μL (1 $\mu\text{L}/\text{min}$ to 10 $\mu\text{L}/\text{min}$)
 - Electrode 10–50 μL (10 $\mu\text{L}/\text{min}$ to 50 $\mu\text{L}/\text{min}$)
- Micro 50–200 μL probe: Flow rate of 50 $\mu\text{L}/\text{min}$ to 200 $\mu\text{L}/\text{min}$. Compatible electrodes are:
 - Electrode 50–200 μL (50 $\mu\text{L}/\text{min}$ to 200 $\mu\text{L}/\text{min}$)
- ESI Calibration probe: Flow rate of 1 $\mu\text{L}/\text{min}$ to 2000 $\mu\text{L}/\text{min}$. Compatible electrodes are:
 - ESI CAL electrode (1 $\mu\text{L}/\text{min}$ to 2000 $\mu\text{L}/\text{min}$)
- Nano < 1 μL probe: Flow rate of 100 nL/min to 1000 nL/min. Compatible electrodes are:
 - Nano electrode (100 nL/min to 1000 nL/min)


Ion Source Components

Micro Probe

Note: Only one Micro probe can be installed at a time. If micro functionality is enabled, then the ESI Calibration probe or a probe port plug must be installed in the front port.

Ion Source Overview


Figure 2-1 Ion Source Components (Micro Probe)

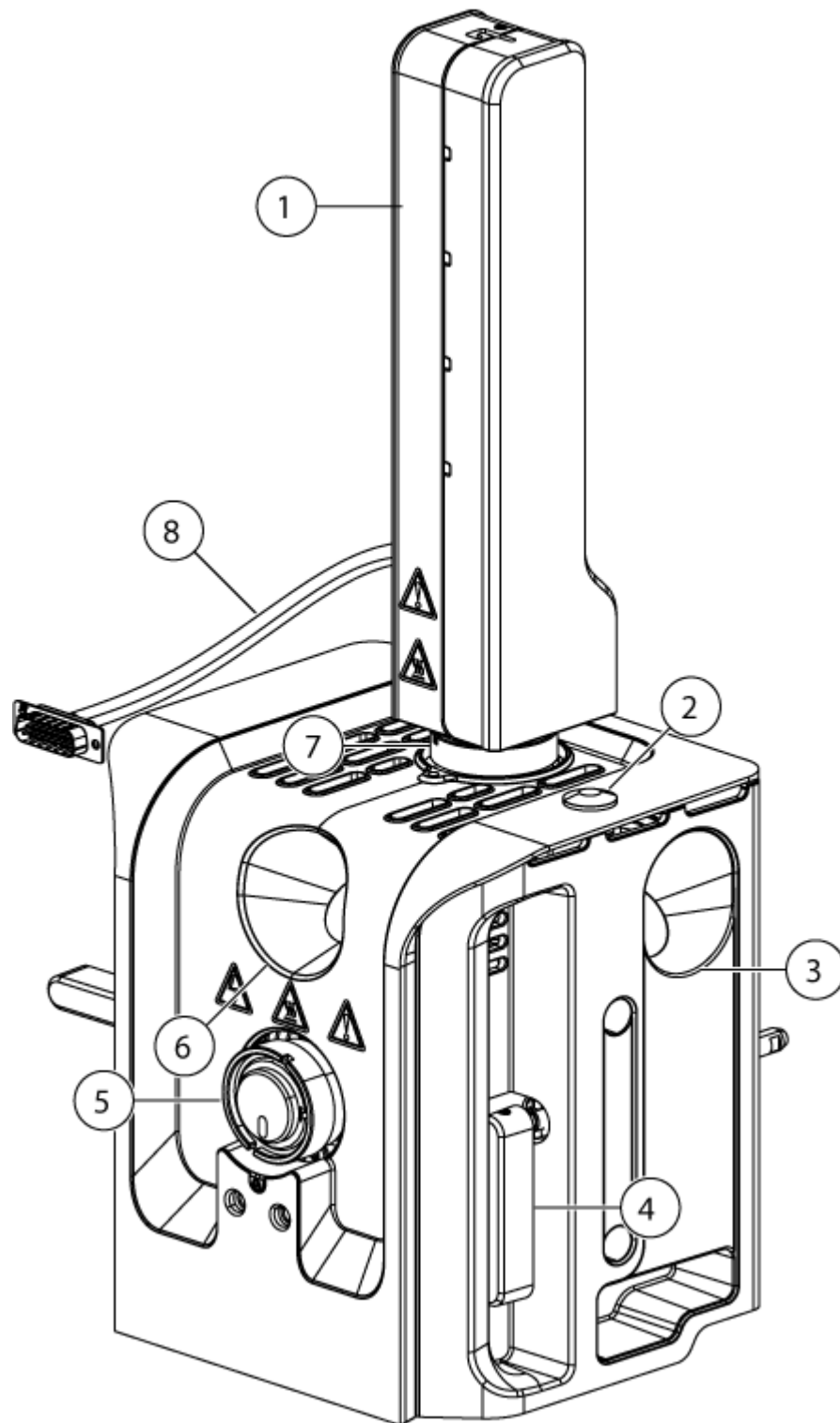

Item	Description
1	Micro column heater. The Micro column heater must be connected to the mass spectrometer, and then configured in the SCIEX OS software. The maximum column heater temperature is 90 °C (194 °F).
2	Top cover plug.
3	Side window port. There is one on each side of the ion source.
4	Source latch. There is one on each side of the ion source. The latches secure the ion source to the mass spectrometer.
5	Front port (ESI Calibration probe or probe port plug). The probe port plug is shown in the figure.
6	Front window port.
7	Top port (Micro probe).
8	Power and communication cable. The cable connects to the SOURCES connection on the mass spectrometer, using an adapter cable included with the ion source.

Nano Probe

Figure 2-2 Ion Source Components (Nano Probe)

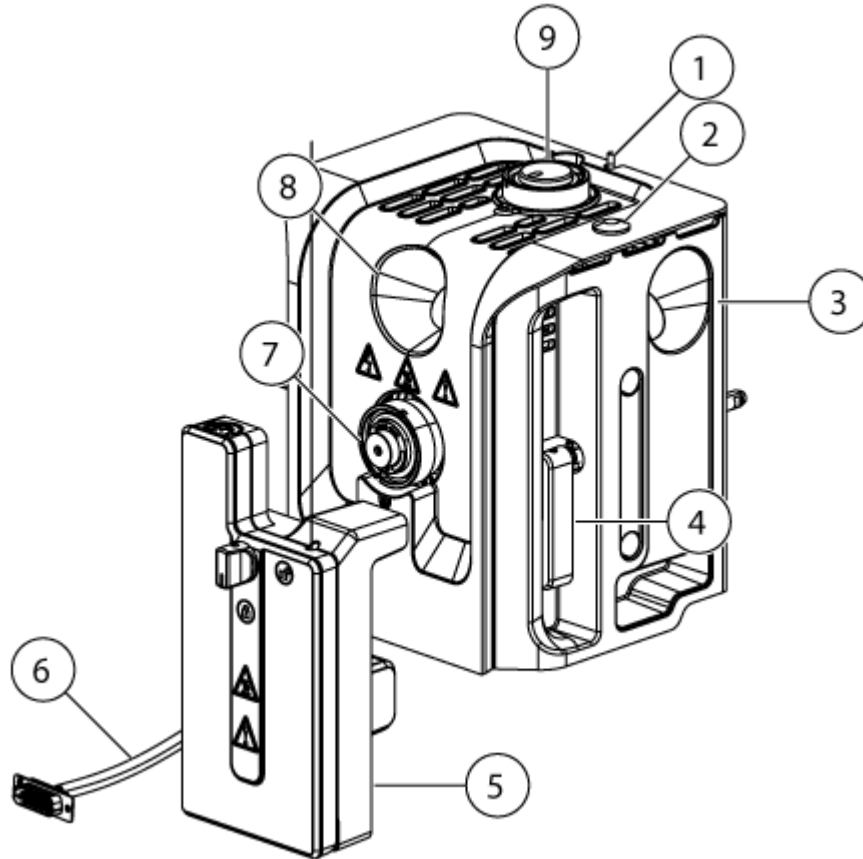

Item	Description
1	High-voltage enable switch. The switch allows the firmware to supply power to the ion source when activated.
2	Top cover plug.
3	Side window port. There one on each side of the ion source.
4	Source latch. There is one on each side of the ion source. The latches secure the ion source to the mass spectrometer.
5	Nano column heater. The Nano column heater must be connected to the mass spectrometer, and then configured in the SCIEX OS software. The maximum column heater temperature is 90 °C (194 °F).


Item	Description
6	Power and communication cable. The cable connects to the SOURCES connection on the mass spectrometer, using an adapter cable included with the ion source.
7	Front port (Nano probe).
8	Front window port.
9	Top port. The probe port plug must be installed for Nano flow operation.

Probe for the Echo[®] MS System

Note: The ZenoTOF 7600 system supports the Echo[®] MS+ system.

Note: Install the probe for the Echo[®] MS system on the top port and the ESI calibration probe in the front port of the OptiFlow Turbo V source.

Figure 2-3 Ion Source Components (Echo[®] MS+ System)

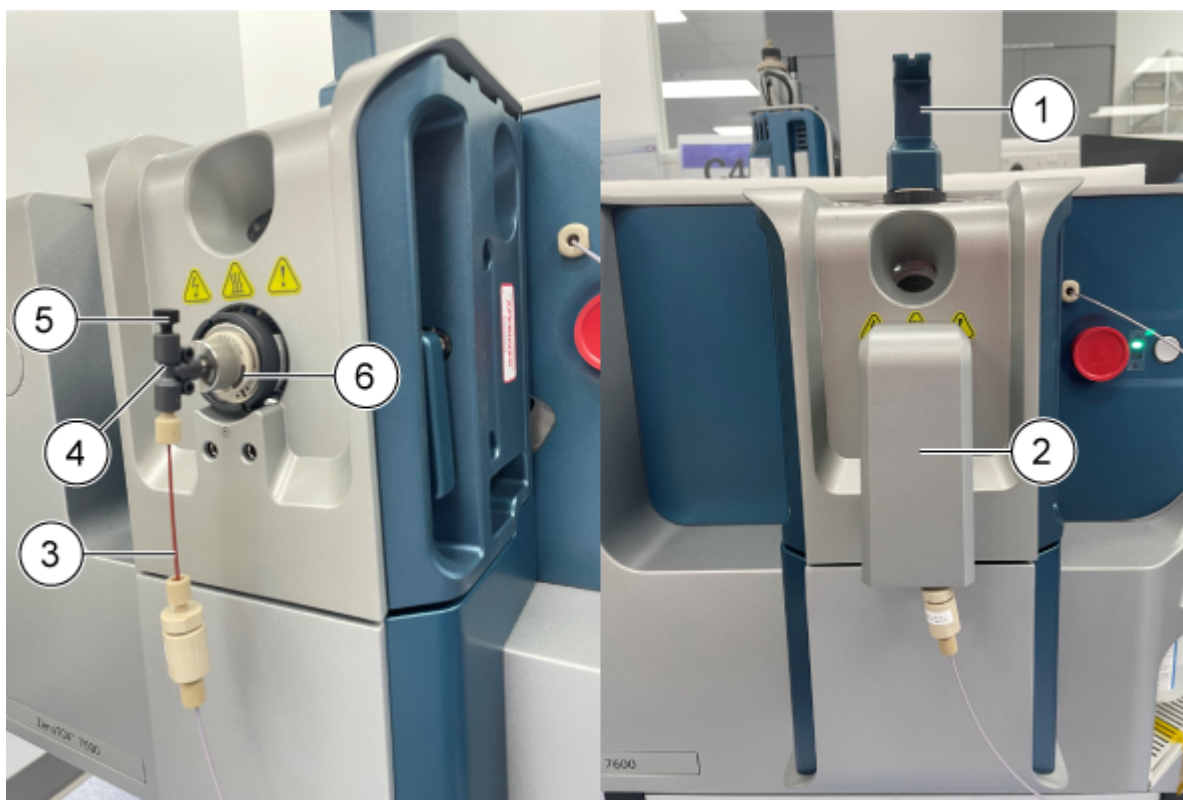

Ion Source Overview

Item	Description
1	Infusion adapter
2	Protective cap
3	CDS tubing
4	T-fitting
5	Plug
6	ESI calibration probe

Gas and Electrical Connections

Gas and low- and high-voltage electrical connections are provided on the front plate of the vacuum interface and they connect internally through the ion source housing. When the ion source is installed on the mass spectrometer, all of the electrical and gas connections are complete.

Ion Source Sense Circuit

An ion source sense circuit disables the high-voltage power supply for the mass spectrometer and the source exhaust system under these conditions:

- The ion source is not installed or is incorrectly installed.
- A probe or plug is not installed or is incorrectly installed.
- The mass spectrometer senses a gas fault.
- The ion source has overheated.
- The infusion adapter is not installed.

Source Exhaust System



WARNING! Ionizing Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Make sure that the source exhaust system is connected and functioning, to safely remove sample vapor exhaust from the laboratory environment. Emissions from the equipment must be exhausted in the general building exhaust and not allowed to exhaust in the workspace of the laboratory. For requirements for the source exhaust system, refer to the document: *Site Planning Guide*.



WARNING! Ionizing Radiation Hazard, Biohazard, or Toxic Chemical Hazard. Vent the source exhaust system to either a dedicated laboratory fume hood or an external ventilation system to prevent hazardous vapors from being released in the laboratory environment.



WARNING! Ionizing Radiation Hazard, Biohazard, or Toxic Chemical Hazard. If an LC system is used with the mass spectrometer, and if the source exhaust system is not functioning properly, then shut down the LC system until functionality of the source exhaust system has been restored.



WARNING! Fire Hazard. Do not direct more than 2 mL/min of flammable solvent in the ion source. Exceeding the maximum flow rate can cause solvent to accumulate in the ion source. Do not use the ion source if the source exhaust system is not enabled and functioning when the ion source and the probe are properly installed.

Note: Make sure that all of the exhaust tubing is securely connected to decrease the chance of equipment exhaust entering the room.

An ion source produces both sample and solvent vapors. These vapors are a potential hazard to the laboratory environment. The source exhaust system is designed to help the user safely remove and correctly handle the sample and solvent vapors. When the ion source is installed, the mass spectrometer does not operate unless the source exhaust system is operating.

An active exhaust system removes ion source exhaust, including gases, solvent, and sample vapor, through a drain port, without introducing chemical noise. The drain port connects through a drain chamber and a source exhaust pump to a drain bottle, and from there to a customer-supplied exhaust ventilation system. For information about the ventilation requirements for the source exhaust system, refer to the mass spectrometer document: *Site Planning Guide*.

Note: Examine the source exhaust system periodically to make sure that the exhaust tubing is intact and that exhaust is not leaking in the room.

Ion Source Installation

3



WARNING! Electrical Shock Hazard. Do not use electrically conductive tubing or fittings, such as stainless steel or any other metal or metallic compound, with the ion source. A static shock or equipment malfunction might occur. Use only non-electrically conductive tubing and fittings, such as PEEK or PEEK-clad fused silica.



WARNING! Electrical Shock Hazard. When introducing a sample through infusion, remove the infusion adapter to disable the high voltage before inspecting any fittings and tubing for leaks. Contact with liquid leaking from the probe fittings or tubing can result in a static shock if there is high voltage present.



WARNING! Electrical Shock Hazard. Install the ion source on the mass spectrometer as the last step in this procedure. High voltage is present when the ion source is installed.

CAUTION: Potential System Damage. Do not lift or carry the ion source with one hand. The ion source is designed to be lifted or carried using the molded grips on each side of the ion source.

The ion source is connected to the vacuum interface and is held in position by two source latches. The interior of the ion source is visible through the tempered glass window on the front of the ion source.

When the ion source is installed, the software recognizes the ion source and shows the ion source identification.

Required Materials

- Ion source
- Red PEEK tubing (0.005 inch bore)
- Micro, Nano, or ESI Calibration probe
- Probe for the Echo[®] MS system
- Column
- Infusion adapter and PEEK Tee
- Electrode
- PEEK-clad fused silica tubing
- Micro components:
 - Micro column heater
 - Top and bottom fittings for the Micro probe
- Nano components:
 - Nano column heater
 - Union and fitting for the Nano electrode

Install the Probe for the Echo[®] MS System or the Micro Probe On the Ion Source

CAUTION: Potential System Damage. Do not let the protruding electrode touch any part of the ion source housing, to avoid damaging the electrode.


1. Slowly insert the micro probe or the probe for the Echo[®] MS system to the top port of the ion source.

Note: If the Nano probe is inserted, then make sure that the nano probe is inserted on the front port of the ion source.

2. Align the dot on the probe with the corresponding dot on the ion source housing.

Ion Source Installation

Figure 3-1 Align the Dot On the Probe

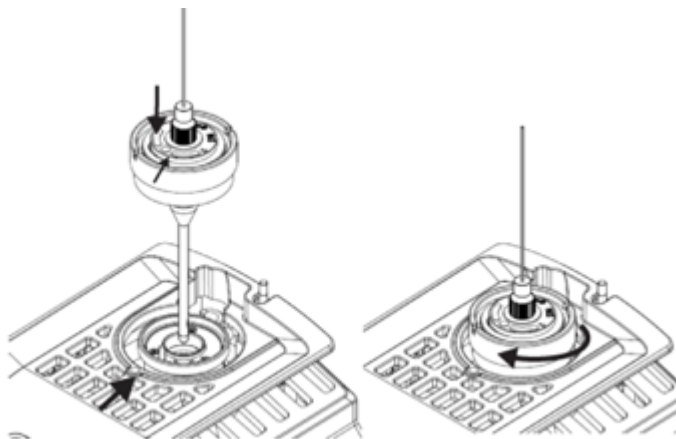

3. Slowly turn the probe in the clock wise direction until it is fully tightened.
4. Install the electrode. Refer to the section: [Connect the Echo[®] MS+ System to the OptiFlow Turbo V Ion Source](#).

Note: Install the front port plug with the probe for the Echo[®] MS system.

Figure 3-2 Front Port Plug



Item	Description
1	Front port plug

Install the Ion Source on the Mass Spectrometer

CAUTION: Potential System Damage. Do not let the protruding electrode touch any part of the ion source housing, to avoid damaging the electrode.

1. Make sure that the source latches on either side of the ion source are pointing up in the 12 o'clock position.
2. Align the ion source with the vacuum interface, making sure that the guide pins on the ion source are aligned with the sockets in the vacuum interface.

Ion Source Installation

3. Push the ion source gently against the vacuum interface, and then rotate the ion source latches down to lock the ion source in position.

Install the Electrode in a Micro Probe



WARNING! Electrical Shock Hazard. Do not use electrically conductive tubing or fittings, such as stainless steel or any other metal or metallic compound, with the ion source. A static shock or equipment malfunction might occur. Use only non-electrically conductive tubing and fittings, such as PEEK or PEEK-clad fused silica.



WARNING! Electrical Shock Hazard. When introducing a sample through infusion, remove the infusion adapter to disable the high voltage before inspecting any fittings and tubing for leaks. Contact with liquid leaking from the probe fittings or tubing can result in a static shock if there is high voltage present.



WARNING! Puncture Hazard. Be careful when handling the electrode. The tip of the electrode is extremely sharp.

CAUTION: Potential System Damage. Install the probe in the ion source before the electrode is installed in the probe. This reduces the risk of damaging the electrode tip while installing it in the ion source.

CAUTION: Potential System Damage. Do not let the protruding electrode tip touch any part of the ion source housing, to avoid damaging the electrode.

Prerequisite Procedures

- [Install the Probe for the Echo® MS System or the Micro Probe On the Ion Source](#)
- [Install the Ion Source on the Mass Spectrometer.](#)

Figure 3-3 Ion Source Probes




Item	Description	Comments
1	Micro probe	The probe is ready to have the electrode installed.
2	ESI Calibration probe or a probe port plug	The probe port plug is displayed in the figure.

1. Install the electrode in the probe, inserting the fused silica or steel end first.
2. Rotate the electrode slightly to install against it in the probe, and then make sure that the tip of the electrode is visible below the end of the probe.
The nominal protrusion for the electrode is 1.0 mm.

Ion Source Installation


Figure 3-4 Electrode in the Micro Probe

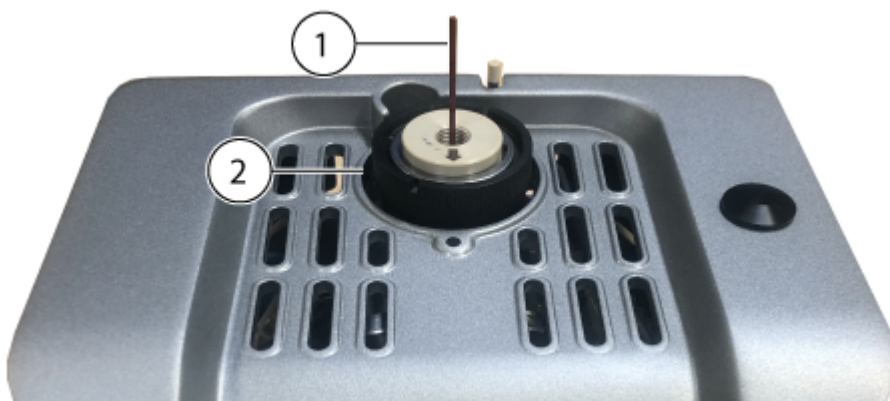

Item	Description
1	Electrode
2	Micro probe

3. Install the bottom fitting over the electrode, and then tighten the fitting until it is finger tight.


Figure 3-5 Bottom Fitting Installed

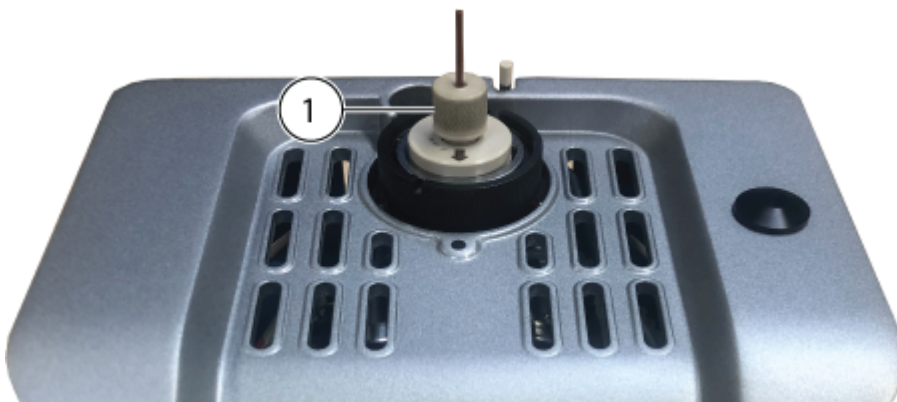

Item	Description
1	Bottom fitting

4. Put the PEEK ferrule on the top fitting, and then put the top fitting on top of the bottom fitting.
The top fitting fits loosely on top of the bottom fitting, to accommodate different column fitting depths.

Figure 3-6 Top Fitting

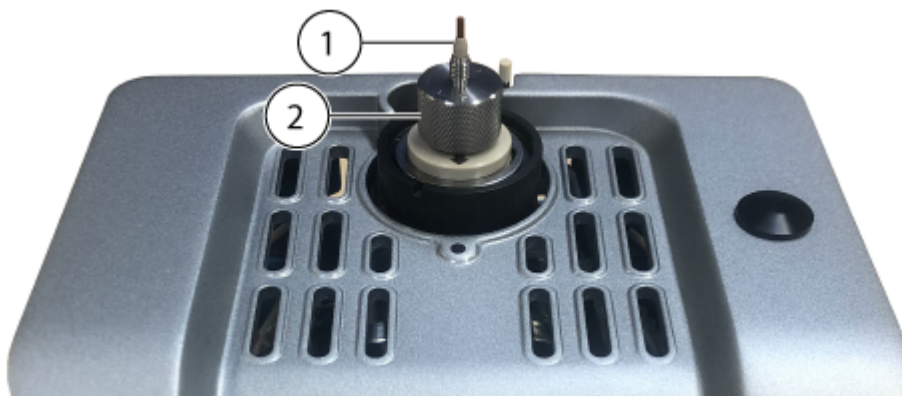

Item	Description
1	PEEK ferrule
2	Top fitting

The electrode installation is complete. The sample tubing, either a column or an infusion adapter and PEEK Tee, can be installed. To install the column, refer to the section: [Install the Micro Column and Heater](#). To install an infusion adapter and PEEK Tee, refer to the section: [Install an Infusion Adapter and Connect an Infusion Line to a Micro Probe](#).

Install an Infusion Adapter and Connect an Infusion Line to a Micro Probe



WARNING! Electrical Shock Hazard. Do not use electrically conductive tubing or fittings, such as stainless steel or any other metal or metallic compound, with the ion source. A static shock or equipment malfunction might occur. Use only non-electrically conductive tubing and fittings, such as PEEK or PEEK-clad fused silica.



WARNING! Electrical Shock Hazard. When introducing a sample through infusion, remove the infusion adapter to disable the high voltage before inspecting any fittings and tubing for leaks. Contact with liquid leaking from the probe fittings or tubing can result in a static shock if there is high voltage present.

Ion Source Installation

Prerequisite Procedures


- [Install the Ion Source on the Mass Spectrometer.](#)
- [Install the Electrode in a Micro Probe.](#)

Samples can be introduced through a direct connection to the electrode to optimize the ion source and mass spectrometer. This is done by using a PEEK union for direct infusion from a syringe pump, or by using a PEEK Tee to combine syringe pump flow with LC mobile phases, such as Tee infusion. Tee infusion is effective for ion source optimization because the solvent composition can be adjusted through the LC system to be similar to the LC elution composition of the subject analyte, thus achieving optimal system performance.

Note: This procedure describes Tee infusion. For direct infusion, substitute a PEEK union for the PEEK Tee.

1. Put the top fitting on top of the bottom fitting. Install the top fitting in the PEEK Tee, and then hold the PEEK Tee and the top fitting clockwise until the fitting is finger tight.


Figure 3-7 PEEK Tee

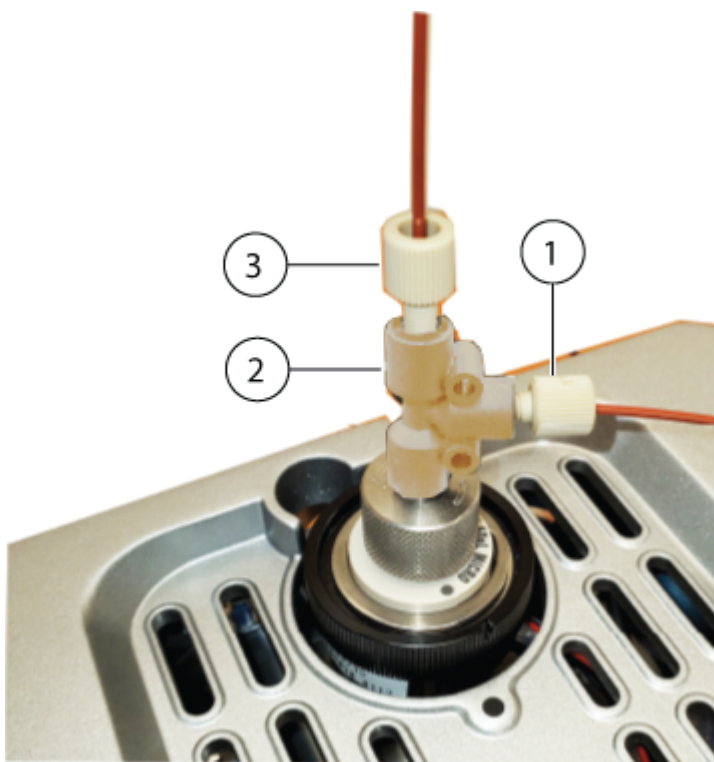

Item	Description
1	Infusion inlet
2	PEEK Tee
3	Mobile phase inlet from the LC system

2. Connect the mobile phase line from the LC system to one inlet of the Tee.
3. Install the infusion line between the infusion inlet and the syringe pump.
4. To make the system operational, install the infusion adapter on the ion source. Install the adapter mounting post in the position hole on the ion source to enable operation of the system. Refer to the figure: [Figure 3-10](#).

Figure 3-8 Infusion Adapter

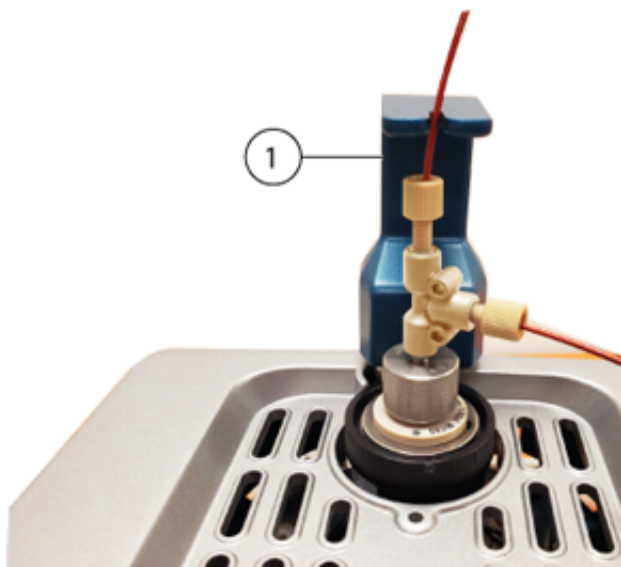

Item	Description
1	Infusion adapter

Install the Micro Column and Heater



WARNING! Electrical Shock Hazard. Make sure that the ion source is completely disconnected from the mass spectrometer before proceeding.

Ion Source Installation



WARNING! Hot Surface Hazard. Beware of burns. The column can become hot during operation. Allow the column to cool before removing it or replacing the PEEK clad fused silica tubing.



WARNING! Electrical Shock Hazard. Do not use electrically conductive tubing or fittings, such as stainless steel or any other metal or metallic compound, with the ion source. A static shock or equipment malfunction might occur. Use only non-electrically conductive tubing and fittings, such as PEEK or PEEK-clad fused silica.




WARNING! Electrical Shock Hazard. When introducing a sample through infusion, remove the infusion adapter to disable the high voltage before inspecting any fittings and tubing for leaks. Contact with liquid leaking from the probe fittings or tubing can result in a static shock if there is high voltage present.

Prerequisite Procedures

- [Install the Ion Source on the Mass Spectrometer.](#)
- [Install the Electrode in a Micro Probe.](#)

1. Install the sample tubing between the column and the LC system. Use the sample tubing provided with the LC system. Refer to the LC System document: *Operator Guide*.
2. Attach the column to the top fitting on the probe, and then tighten the fitting until it is finger tight. Make sure that the electrode is fully seated in the column fitting to minimize the possibility of any dead volume. Hold the column, and then rotate the top fitting counterclockwise until it is finger tight.

Figure 3-9 Column

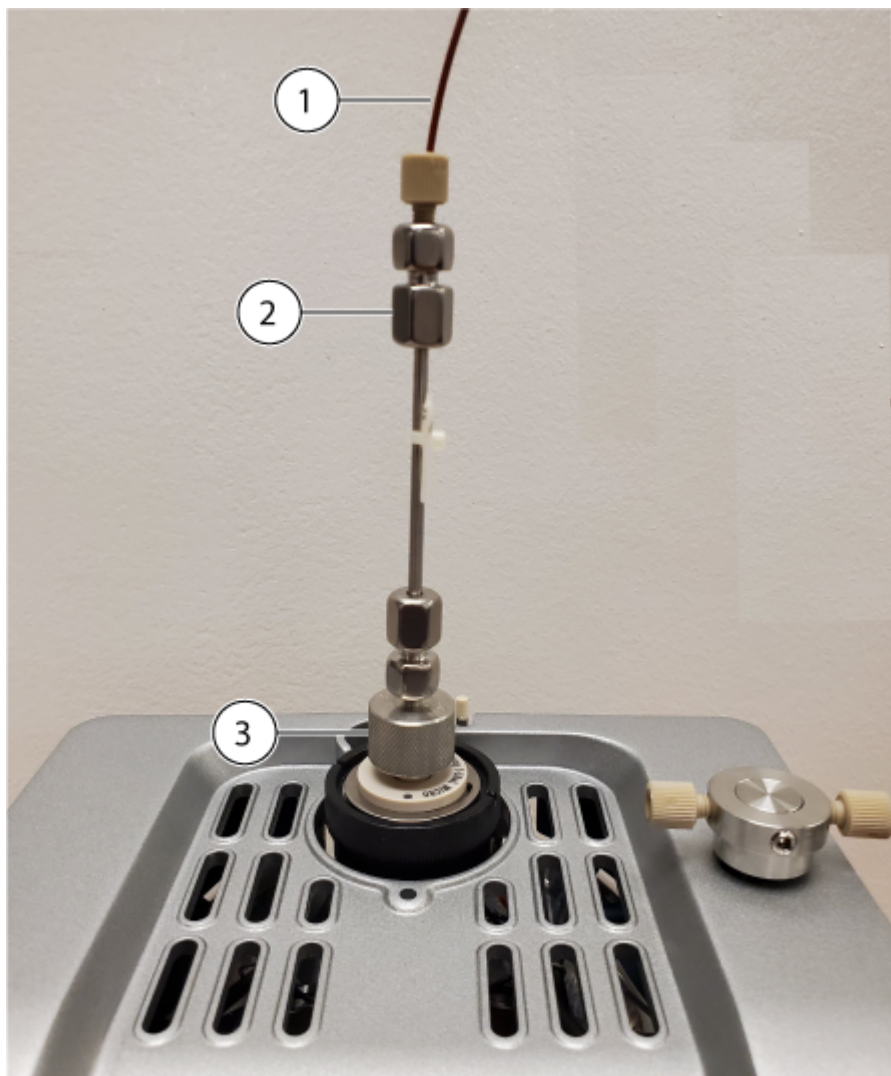


Item	Description
1	Sample tubing
2	Column
3	Top fitting

3. Install the mounting post of the column heater in the position hole on the ion source.

Ion Source Installation


Figure 3-10 Position Hole

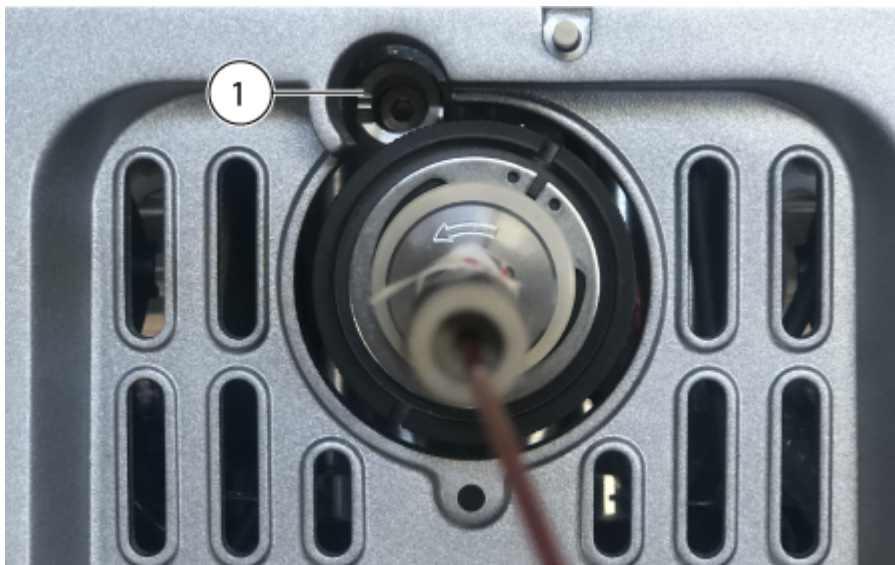

Item	Description
1	Position hole for the column heater mounting post

4. Rotate the left side of the column heater toward the column.

Figure 3-11 Column Heater, Left Side

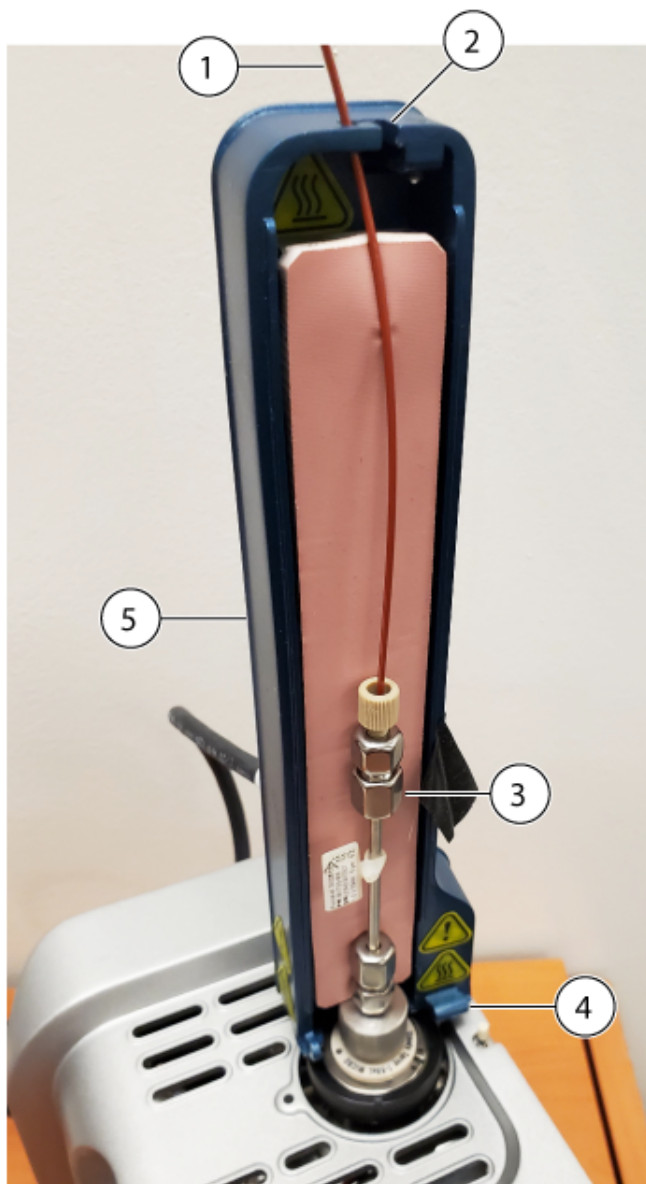

Item	Description
1	Sample introduction tubing
2	Guide slot for the sample tubing
3	Column
4	Hinge


Ion Source Installation

Item	Description
5	Left side of column heater Note: The heater is made of two parts which must be assembled around the column.

Make sure that the mounting post is seated firmly in the position hole on the ion source.

5. Route the PEEK-clad fused silica tubing through the entrance slot at the top of the column heater. Refer to the figure: [Figure 3-11](#).
6. Put the right side of the column heater on the hinge at the base of the left side of the column heater, and then close both sides of the heater until they lock together.

Figure 3-12 Column Heater

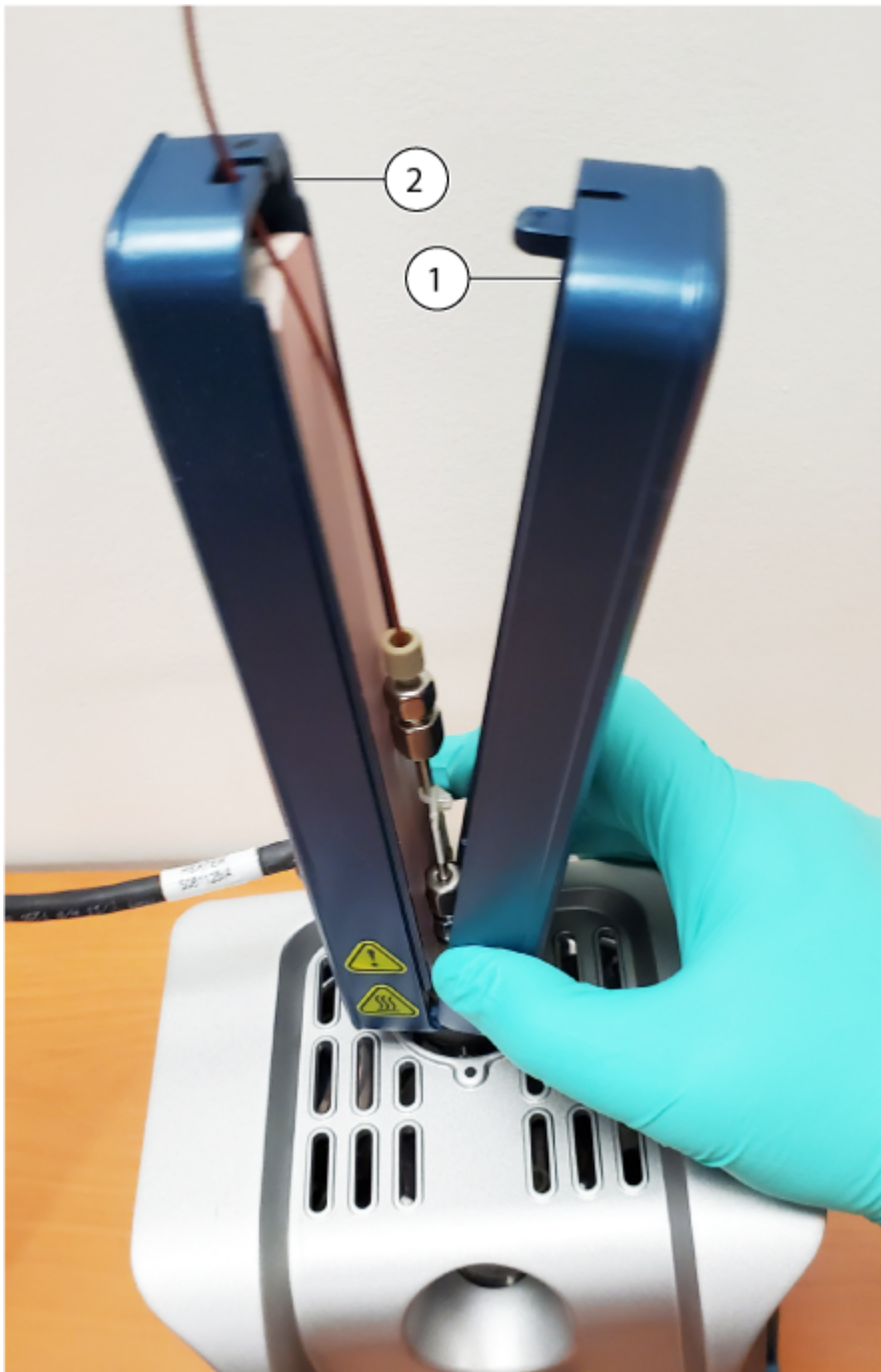

Item	Description
1	Right side of the column heater

Ion Source Installation

Item	Description
2	Left side of the column heater

7. Connect the column heater power and communication cable to the **SOURCES** connection on the mass spectrometer, using an adapter cable included with the ion source.
8. Configure the column heater in the SCIEX OS software.

Install the Electrode in the Nano Probe



WARNING! Electrical Shock Hazard. Do not use electrically conductive tubing or fittings, such as stainless steel or any other metal or metallic compound, with the ion source. A static shock or equipment malfunction might occur. Use only non-electrically conductive tubing and fittings, such as PEEK or PEEK-clad fused silica.



WARNING! Puncture Hazard. Be careful when handling the electrode. The tip of the electrode is extremely sharp.

CAUTION: Potential System Damage. Install the probe in the ion source before the electrode is installed in the probe. This reduces the risk of damaging the electrode tip while installing it in the ion source.

CAUTION: Potential System Damage. Do not let the protruding electrode tip touch any part of the ion source housing, to avoid damaging the electrode.

Prerequisite Procedures

- [Install the Ion Source on the Mass Spectrometer.](#)

1. Install the fitting in the union and then rotate the fitting clockwise until it is finger tight.

Figure 3-13 Nano Electrode



Item	Description	Comments
1	Union	—
2	Fitting	The fitting is threaded into the union
3	Nano electrode	—

Note: The fitting is pre-installed on the electrode.

2. If the front probe port has a probe port plug installed, then remove the plug and then install the Nano probe.
 - a. Loosen the knurled ring on the probe port plug, and then gently pull the plug straight out of the ion source.
 - b. Install the Nano probe in the front probe port.

Tip! When installing the probe in the ion source, align the dot on the probe with the corresponding dot on the ion source housing.

- c. Tighten the knurled ring on the Nano probe.

Figure 3-14 Nano Probe



Item	Description
1	Nano probe
2	Dot on the probe Note: The dot on the probe shows the location of a pin on the ion source housing that fits in a hole in the back of the probe. When the dot on the probe is aligned with the dot on the ion source housing, the pin and the hole are aligned correctly.
3	Dot on ion source housing

3. Install the union, with the electrode attached, in the Nano probe, and then rotate the union clockwise to tighten it.

Figure 3-15 Union Installed in the Nano Probe



Item	Description
1	Nano probe
2	Union

Install the Nano Column Cartridge and Heater



WARNING! Electrical Shock Hazard. Make sure that the ion source is completely disconnected from the mass spectrometer before proceeding.



WARNING! Hot Surface Hazard. Beware of burns. The column can become hot during operation. Allow the column to cool before removing it or replacing the PEEK clad fused silica tubing.



WARNING! Electrical Shock Hazard. Do not use electrically conductive tubing or fittings, such as stainless steel or any other metal or metallic compound, with the ion source. A static shock or equipment malfunction might occur. Use only non-electrically conductive tubing and fittings, such as PEEK or PEEK-clad fused silica.

Ion Source Installation



WARNING! Electrical Shock Hazard. Do not connect the Nano column heater to a European style two-pin electrical outlet. Fire or electric shock might occur.




WARNING! Personal Injury Hazard. Always wear protective eyewear when handling the Nano column heater. A part of the fused silica tubing protrudes from the top of the column heater and can cause injury.

Prerequisite Procedures

- [Install the Ion Source on the Mass Spectrometer.](#)
- [Install the Electrode in the Nano Probe.](#)

1. If the ion source is to be connected to an LC system, then attach a fitting to the pre-column transfer tube, and then tighten the fitting clockwise in the Nano column cartridge until it is finger-tight.


Figure 3-16 Nano Column Cartridge

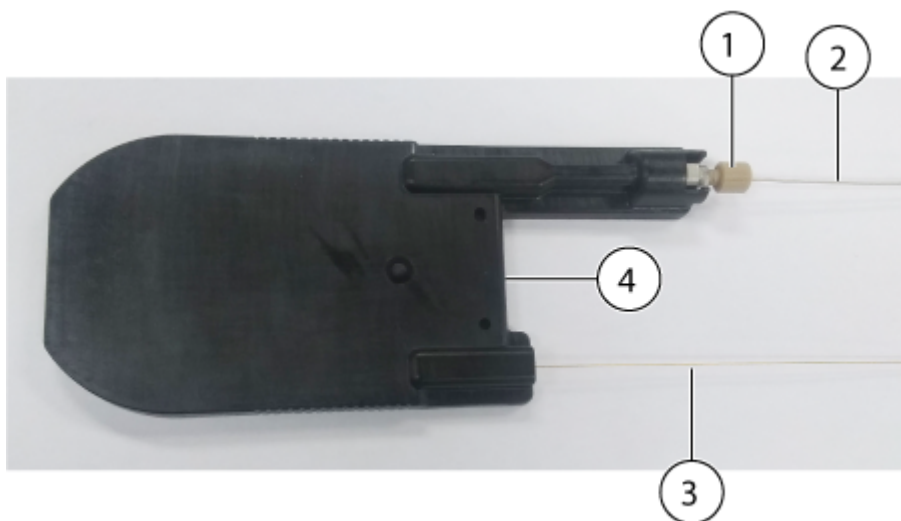

Item	Description
1	Fitting
2	Pre-column transfer tube, connected to an LC system
3	Post-column transfer tube, connected to the Nano probe
4	Nano column cartridge

2. Position the column cartridge in the back part of the Nano column heater.

Figure 3-17 Back Part of the Nano Column Heater

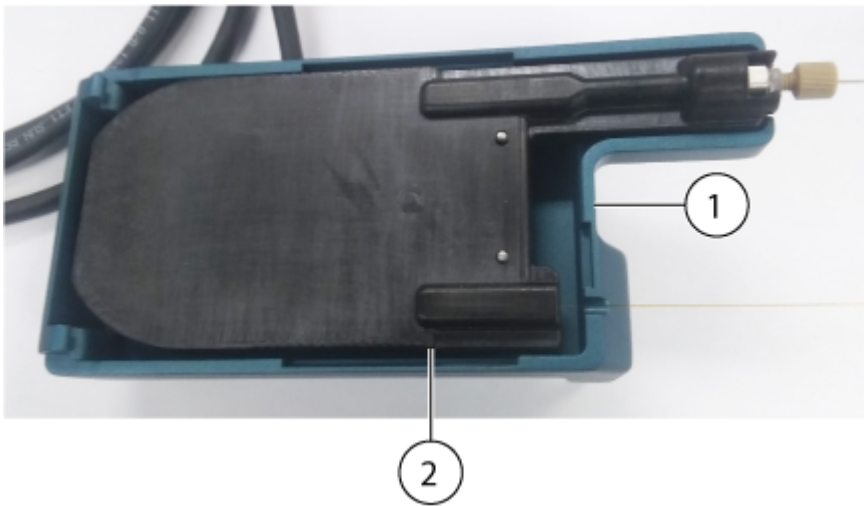

Item	Description
1	Nano column heater
2	Column cartridge

3. Position the base of the front part of the column cartridge over the hinges of the back part, and then fold the front part down to close the column heater.

Ion Source Installation


Figure 3-18 Front and Back Parts of the Nano Column Heater



Item	Description
1	Front part of the Nano column heater
2	Back part of the Nano column heater
3	Locking dial in unlocked position

4. Rotate the locking dial clockwise to secure the column cartridge in the Nano column heater.

Figure 3-19 Nano Column Heater Closed

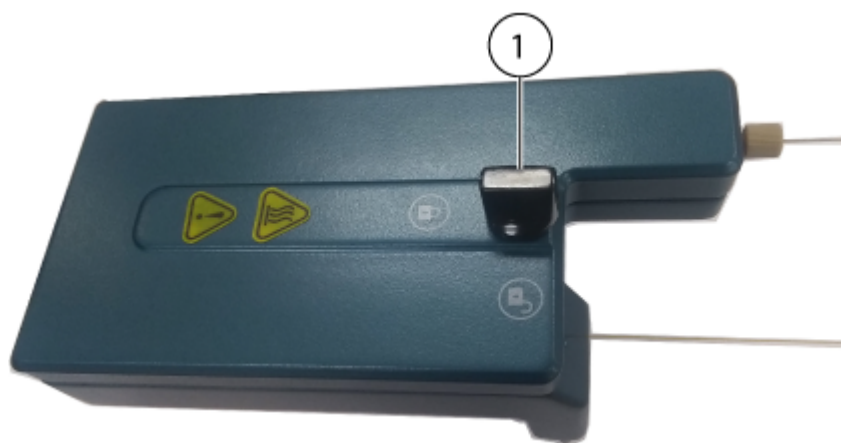

Item	Description
1	Locking dial in locked position

- Find the connection points on the front of the ion source, and then install the Nano column heater on the ion source. Gently press the column heater onto the ion source. Do not use excessive force.

Note: The column heater has two pins that fit in the connection points.

Figure 3-20 Nano Column Heater Connection on the Ion Source




Item	Description
1	Nano probe
2	Connection points for the Nano column heater

- Attach a fitting to the post-column transfer tube, install the fitting in the union installed in the Nano probe, and then rotate the fitting clockwise until it is finger tight.

Ion Source Installation

Figure 3-21 Nano Column Heater Installed on the Ion Source

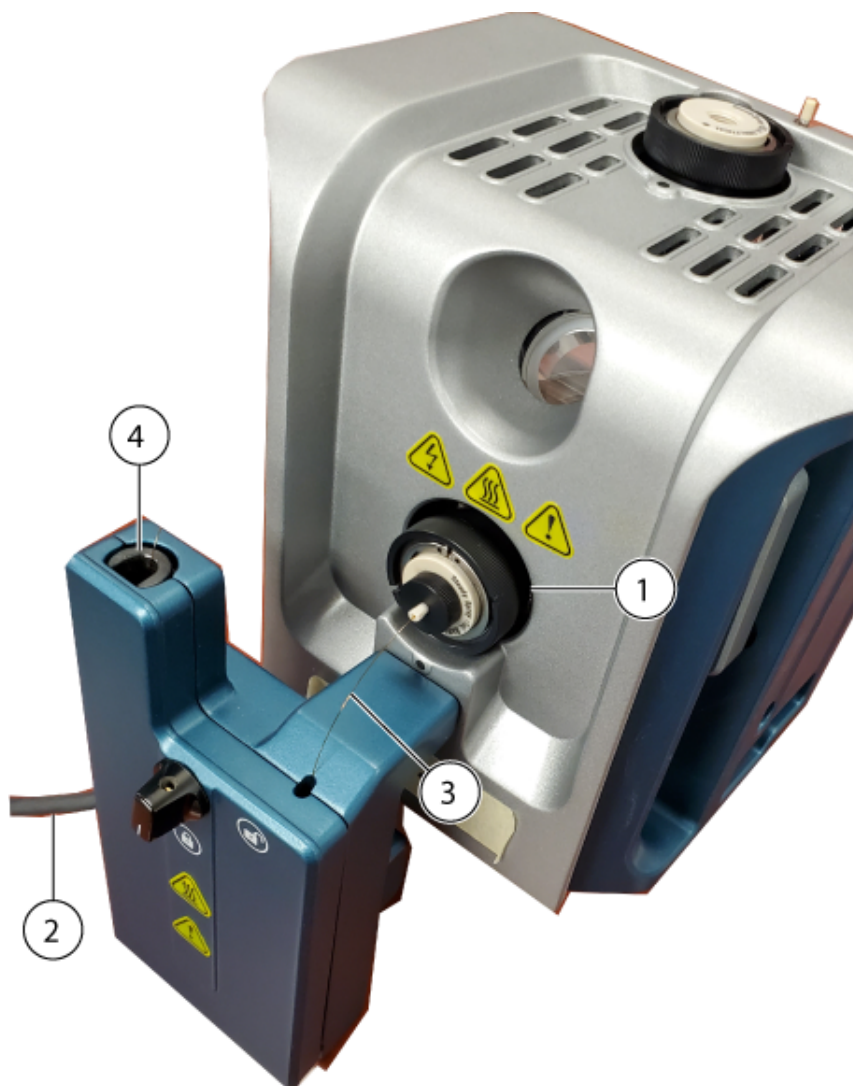

Item	Description
1	Nano probe
2	Power and communication cable to the LC system
3	Post-column transfer tube
4	Pre-column transfer tube port

7. Connect the pre-column transfer tube to the LC system.
8. Connect the column heater power and communication cable to the **SOURCES** connection on the mass spectrometer, using an adapter cable included with the ion source.

-
9. Configure the column heater in the SCIEX OS software.

Connect an Infusion Line to the Nano Probe



WARNING! Electrical Shock Hazard. When introducing a sample through infusion, remove the infusion adapter to disable the high voltage before inspecting any fittings and tubing for leaks. Contact with liquid leaking from the probe fittings or tubing can result in a static shock if there is high voltage present.

Prerequisite Procedures
<ul style="list-style-type: none">• Install the Electrode in the Nano Probe.


Use direct infusion from a syringe pump to introduce a sample. A direct connection to the electrode optimizes the source and mass spectrometer conditions.

- Attach a fitting to the infusion line, and then rotate the fitting clockwise until it is finger tight in the union.

Note: The electrode is installed in the union and then the union is installed in the probe.

Ion Source Installation

Figure 3-22 Infusion Line

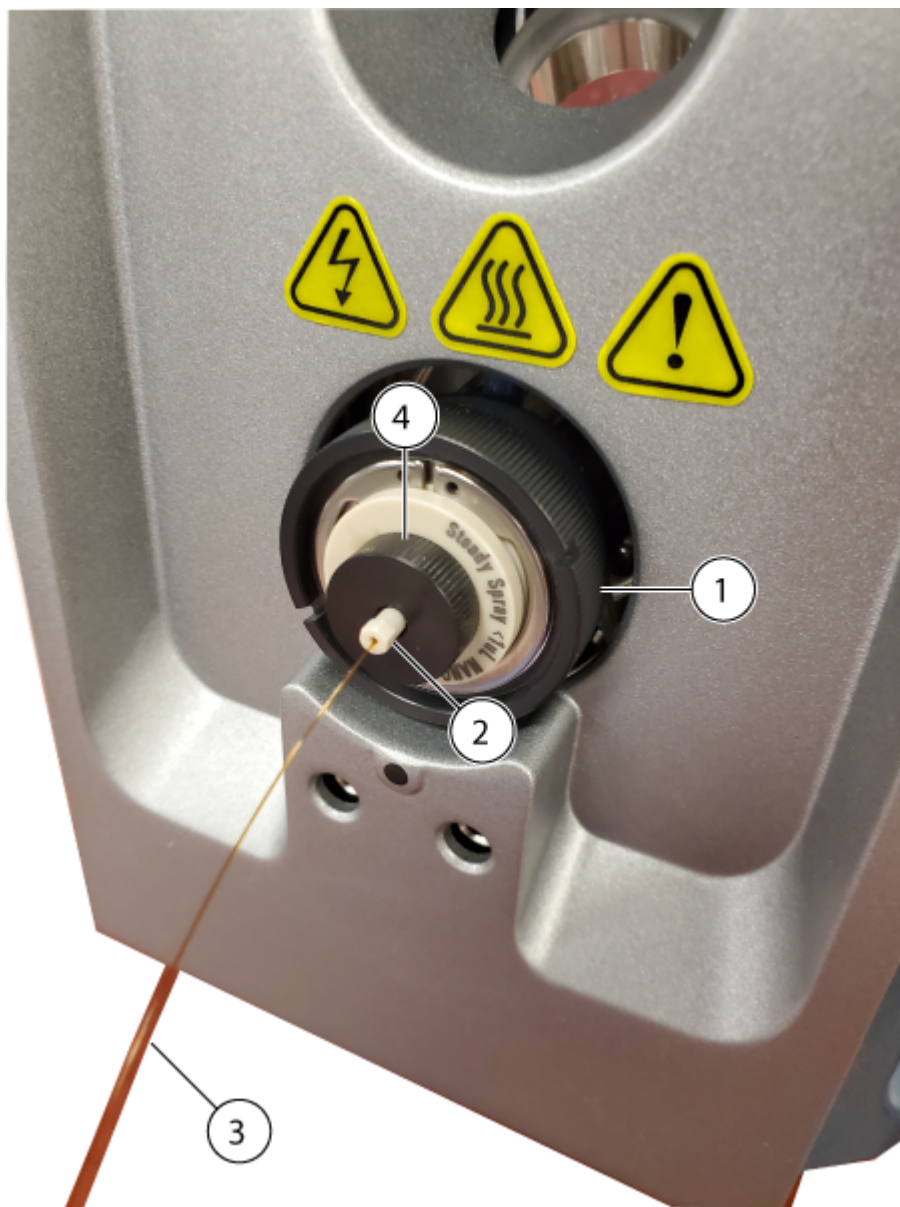


Item	Description
1	Nano probe
2	Fitting
3	Infusion line
4	Union

Install ESI Calibration Probe and Electrode

CAUTION: Potential System Damage. Do not let the protruding electrode touch any part of the ion source housing, to avoid damaging the electrode.

Figure 3-23 Ion Source

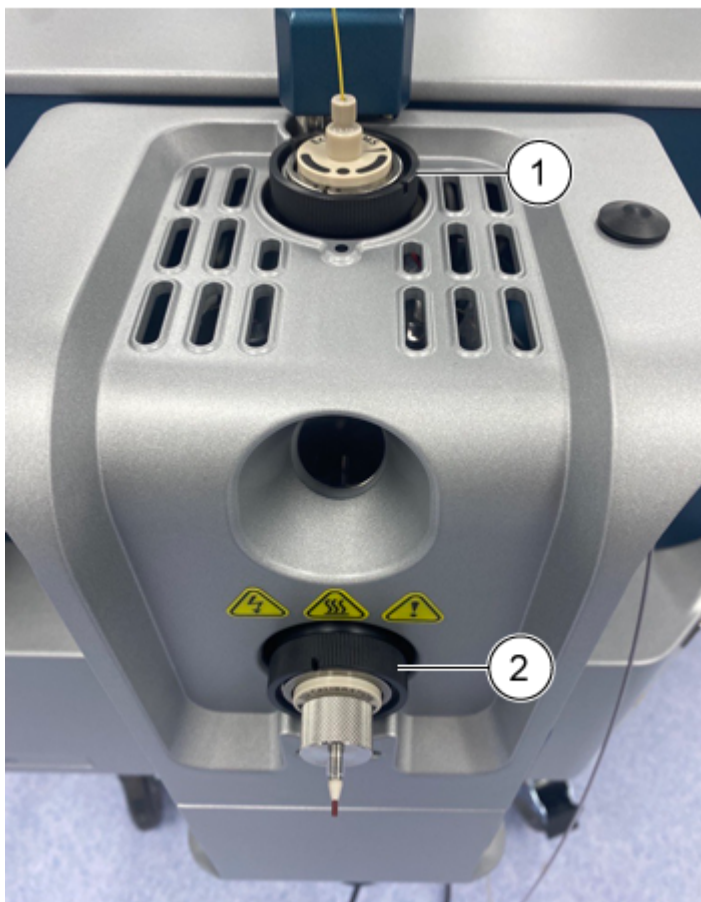


Item	Description
1	Micro probe or the probe for the Echo [®] MS system
2	ESI calibration probe

1. Install the electrode in the ESI calibration probe.

Ion Source Installation

Figure 3-24 ESI Calibration Probe

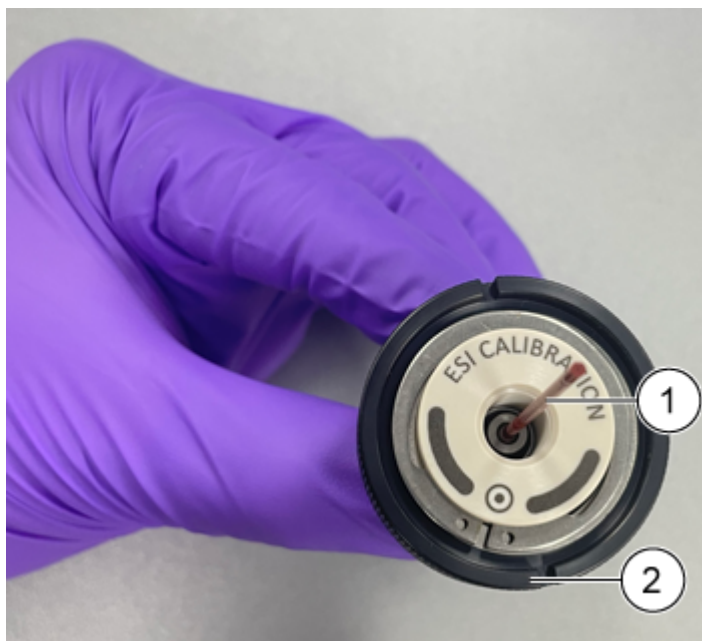


Item	Description
1	Electrode
2	ESI calibration probe

Note: Make sure to install the fused silica or steel end first.

2. Rotate the electrode slightly to install against it in the probe, and then make sure that the tip of the electrode is visible below the end of the probe.
3. Install the bottom fitting over the electrode, and then tighten the fitting until it is finger tight.

Figure 3-25 Bottom Fitting Over the Electrode

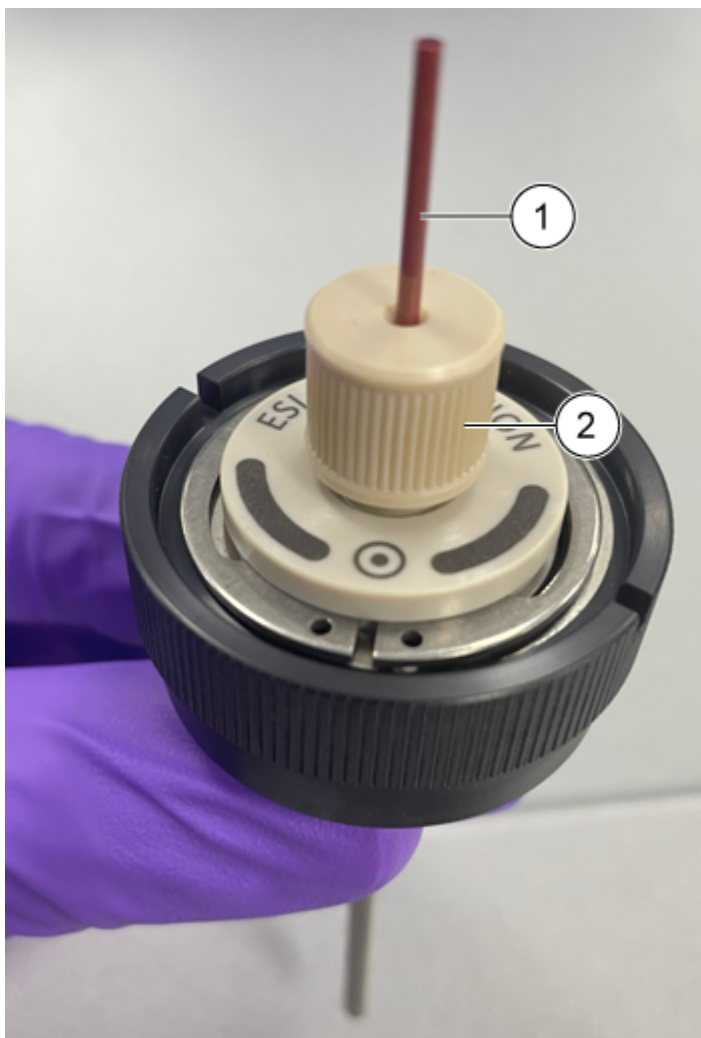

Item	Description
1	Electrode
2	Bottom fitting

- Put the PEEK ferrule on the top fitting, and then put the top fitting on top of the bottom fitting.

Note: To accommodate different column fitting depths, the top fitting fits loosely on top of the bottom fitting.

- Install the ESI calibration probe on the front port of the ion source.

Figure 3-26 ESI Calibration Probe



Item	Description
1	ESI calibration probe
2	Top fitting
3	PEEK Ferrule
4	ESI calibration probe electrode

6. Connect the CDS tubing to the calibration electrode.

Install the Calibration Tubing for the Echo[®] MS+ System



WARNING! Electrical Shock Hazard. Do not use electrically conductive tubing or fittings, such as stainless steel or any other metal or metallic compound, with the ion source. A static shock or equipment malfunction might occur. Use only non-electrically conductive tubing and fittings, such as PEEK or PEEK-clad fused silica.



WARNING! Electrical Shock Hazard. When introducing a sample through infusion, remove the infusion adapter to disable the high voltage before inspecting any fittings and tubing for leaks. Contact with liquid leaking from the probe fittings or tubing can result in a static shock if there is high voltage present.



WARNING! Puncture Hazard. Be careful when handling the electrode. The tip of the electrode is extremely sharp.

CAUTION: Potential System Damage. Do not let the protruding electrode tip touch any part of the ion source housing, to avoid damaging the electrode.


Note: The ZenoTOF 7600+ system does not support the Echo[®] MS+ system.

Prerequisite Procedures

- [Connect the Echo[®] MS+ System to the OptiFlow Turbo V Ion Source](#)
- [Install ESI Calibration Probe and Electrode](#)
- [Install the Ion Source on the Mass Spectrometer](#)

1. Install the plug in the T-fitting, and then tighten it, as shown in the following figure.


Figure 3-27 Assemble the T-fitting

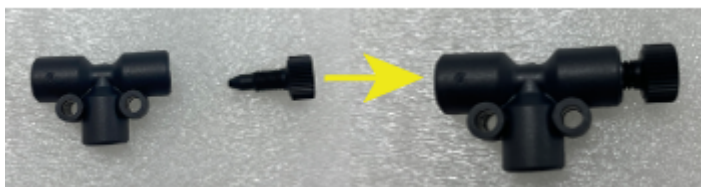

2. Install the T-fitting in the front port of the ion source.

Ion Source Installation

Figure 3-28 T-fitting and CDS Tubing

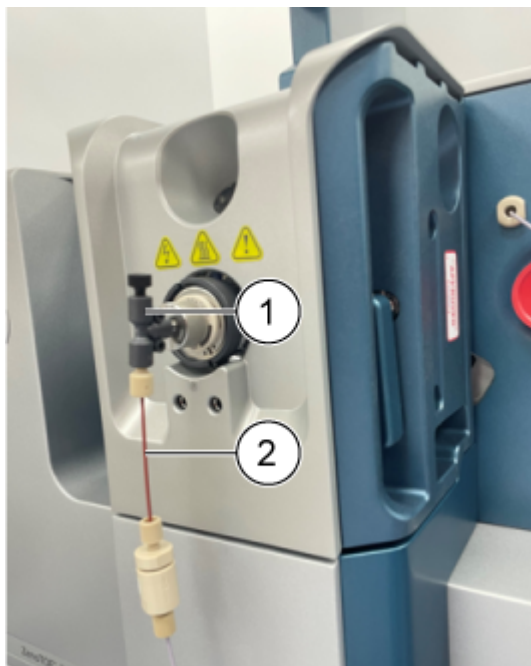

3. Connect the CDS calibration tubing to the bottom of the T-fitting.

Item	Description
1	T-fitting
2	CDS calibration tubing

Note: Make sure that the CDS tubing and the T-fitting are tightly connected.

4. Install the protective cap on the T-fitting.

Figure 3-29 Protective Cap On the T-Fitting



5. To make sure that the system can operate, put the adapter mounting post in the position hole on the ion source.

Figure 3-30 Infusion Adapter



Sample Inlet Requirements

Note: This section is not applicable for the Echo[®] MS+ system. Use applicable analytical procedures and practices to minimize external dead volumes.

- Prefilter samples so that the capillary tubing in the sample inlets is not blocked by particles, precipitated samples, or salts.
- Make sure that all of the connections are tight enough to prevent leaks. Do not over-tighten.

Examine for Leaks



WARNING! Toxic Chemical Hazard. Wear personal protective equipment (PPE), including a laboratory coat, gloves, and safety glasses, to avoid skin or eye exposure.

Regularly examine the ion source for leaks.

1. Make sure that the installed ion source is fully sealed to the mass spectrometer, with no evidence of leaks.
2. Examine fittings and tubing to make sure that there are no leaks.
3. Open the column heater to examine the column connections.

Ion Source Maintenance

4

The following warnings apply to all of the maintenance procedures in this section.



WARNING! Hot Surface Hazard. Before any maintenance procedures are started, let the temperature of the OptiFlow Turbo V ion source decrease for at least 60 minutes. Some surfaces of the ion source and vacuum interface become hot during operation.



WARNING! Fire and Toxic Chemical Hazard. Keep flammable liquids away from flame and sparks and use them only in vented chemical fume hoods or safety cabinets.



WARNING! Toxic Chemical Hazard. Wear personal protective equipment (PPE), including a laboratory coat, gloves, and safety glasses, to avoid skin or eye exposure.



WARNING! Ionizing Radiation Hazard, Biohazard, or Toxic Chemical Hazard. In the event of a chemical spill, review product safety data sheets for specific instructions. Make sure that the system is in Standby state before cleaning a spill near the ion source. Use applicable personal protective equipment and absorbent wipes to contain the spill and dispose of it following local regulations.



WARNING! Electrical Shock Hazard. Avoid contact with the high voltages applied to the ion source during operation. Put the system in Standby state before adjusting the sample tubing or other equipment near the ion source.

CAUTION: Potential System Damage. Do not lift or carry the ion source with one hand. The ion source is designed to be lifted or carried using the molded grips on each side of the ion source.

This section contains general maintenance procedures for the ion source. To identify how often to clean or do maintenance on the ion source, consider the following:

- Compounds tested
- Cleanliness of the samples and sample preparation techniques
- Amount of time an idle probe contains a sample
- Overall system run time

These factors can cause changes in ion source performance, indicating that maintenance is required.

Make sure that the installed ion source is fully sealed to the mass spectrometer with no evidence of gas leaks. Regularly examine the ion source and its fittings for leaks. Clean the ion source components regularly to keep the ion source in good working condition.

Recommended Maintenance Schedule

The following table supplies a recommended schedule for cleaning and maintenance of the ion source. For a list of consumable and spare parts, refer to the document: *Parts and Equipment Guide*.

Tip! Do maintenance tasks regularly to make sure that system performance is optimal.

Contact a Qualified Maintenance Person (QMP) to order consumable parts and for basic service and maintenance requirements. Contact a SCIEX Field Service Employee (FSE) for all other service and maintenance requirements.

Note: For part numbers, refer to the document: *Parts and Equipment Guide*.

Table 4-1 Maintenance Tasks

Component	Frequency	Task	For more information...
Electrode	As needed	Examine and replace	Refer to the section: Install the Electrode in a Micro Probe , or Install the Electrode in the Nano Probe . To replace the electrode for the Echo [®] MS system, refer to the documents: <i>Echo[®] MS+ System User Guide</i> . or the <i>Echo[®] MS System Addendum</i> .
Electrode	As needed	Clean	Refer to the section: Clean the Electrode .
Probe for the Echo [®] MS system	As needed	Examine	Refer to the documents: <i>Echo[®] MS+ System User Guide</i> or the <i>Echo[®] MS System Addendum</i> .
Micro, Nano, or ESI Calibration probes	As needed	Replace	Refer to the section: Remove the Probe .

Ion Source Maintenance

Table 4-1 Maintenance Tasks (continued)

Component	Frequency	Task	For more information...
Sample tubing	As needed	Replace	Refer to the section: Install the Micro Column and Heater or Install an Infusion Adapter and Connect an Infusion Line to a Micro Probe .
Ion source surfaces	As needed	Clean	Refer to the section: Clean the Ion Source Surfaces .


OptiFlow Turbo V Ion Source Handling



WARNING! Personal Injury Hazard. Always wear protective eyewear when handling the Nano column heater. A part of the fused silica tubing protrudes from the top of the column heater and can cause injury.

Surfaces of the ion source become hot during operation. The following figures show surfaces that are cooler (blue) and surfaces that remain hot for an extended period of time (red). Do not touch the surfaces shown in red while using or removing the ion source.

Figure 4-1 OptiFlow Turbo V Ion Source Micro Hot Surfaces (Red=Hot, Blue=Handle with Care)

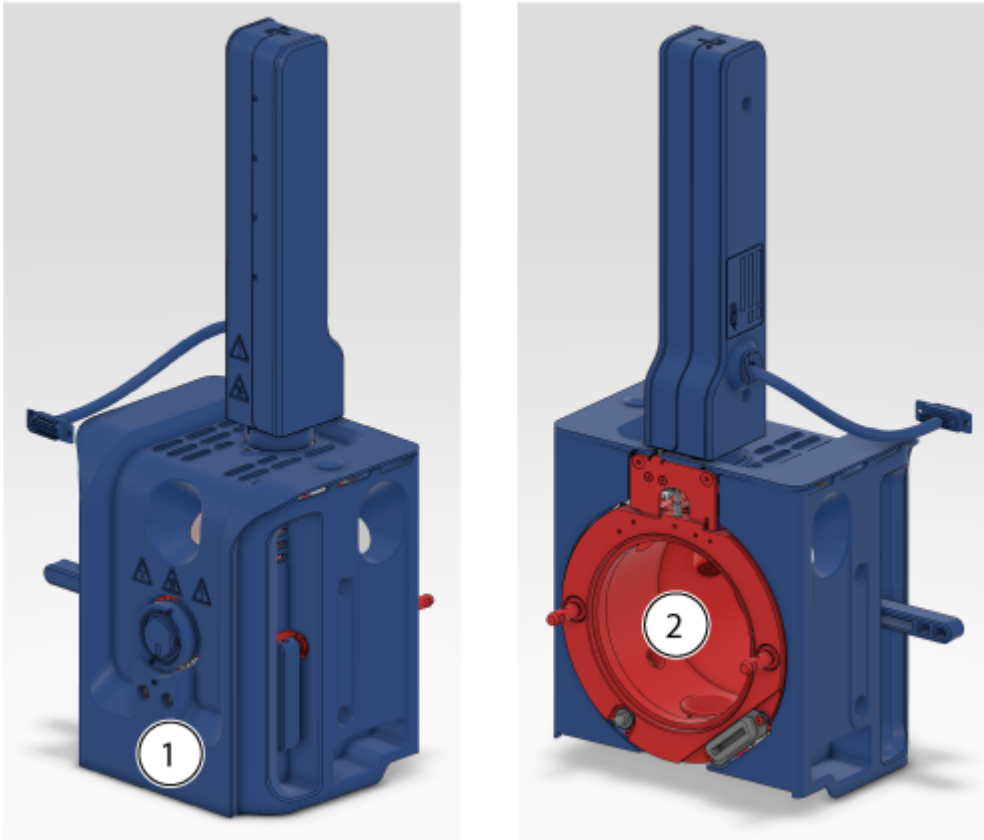

Ion Source Maintenance

Figure 4-2 OptiFlow Turbo V Ion Source Nano Hot Surfaces (Red=Hot, Blue=Handle with Care)


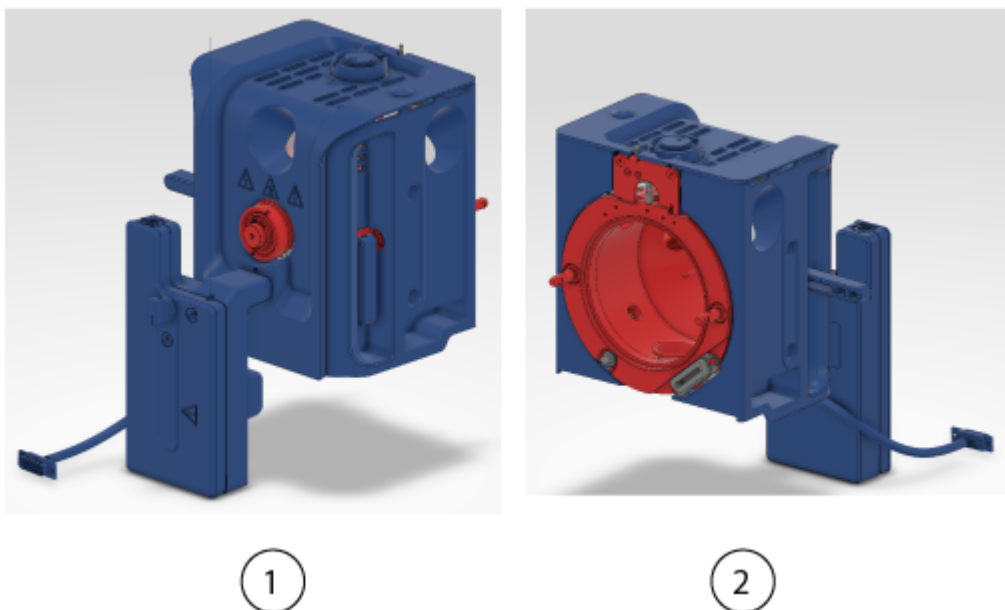


Figure 4-3 OptiFlow Turbo V Ion Source Hot Surfaces for the Echo[®] MS+ System (Red=Hot, Blue=Handle with Care)

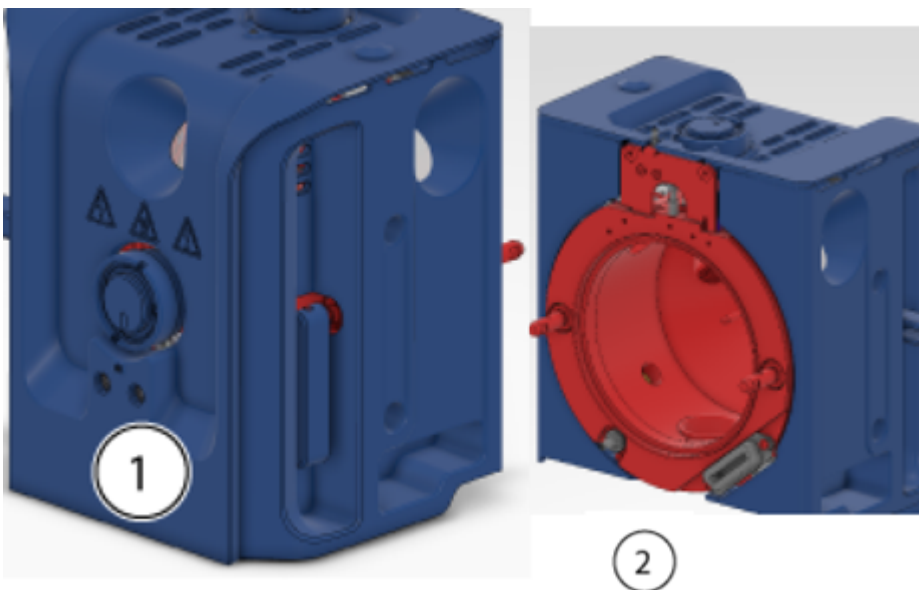

Item	Description
1	Front

Item	Description
2	Back

Remove the Ion Source



WARNING! Hot Surface Hazard. Before any maintenance procedures are started, let the temperature of the OptiFlow Turbo V ion source decrease for at least 60 minutes. Some surfaces of the ion source and vacuum interface become hot during operation.

CAUTION: Potential System Damage. Do not let the protruding electrode tip or the corona discharge needle touch any part of the ion source housing, to avoid damaging the probe.

The ion source can be removed quickly and easily, without tools. SCIEX recommends that the ion source be removed from the mass spectrometer before any maintenance activities are performed.

1. Stop any ongoing scans.
2. Put the mass spectrometer in Standby state.
3. Wait at least 60 minutes for the ion source to cool.
4. If the ion source is configured for use with a column heater, then disconnect the power and communication cable from the **SOURCES** connection on the mass spectrometer.
5. Remove the column and column heater.
 - If the ion source has a Micro column installed, then remove the column heater and disconnect the column from the probe fitting. Refer to the section: [Install the Micro Column and Heater](#).
 - If the ion source has a Nano column cartridge and heater installed, then remove the column heater and disconnect the post-column transfer tube. Refer to the section: [Install the Nano Column Cartridge and Heater](#).
 - If the ion source has the probe for the Echo[®] MS+ system installed, then disconnect the electrode from the Echo[®] MS+ system. Refer to the document: *Echo[®] MS System User Guide*.
6. Disconnect the infusion line.
 - If the ion source has an infusion adapter and PEEK Tee connected to the Micro probe, then disconnect the infusion adapter and PEEK Tee from the probe fitting. Refer to the section: [Install an Infusion Adapter and Connect an Infusion Line to a Micro Probe](#).

Ion Source Maintenance

- If the ion source has an infusion line connected to the Nano probe, then disconnect the infusion line from the probe. Refer to the section: [Connect an Infusion Line to the Nano Probe](#).
7. Rotate the two source latches upward to the 12 o'clock position to release the ion source.
 8. Pull the ion source gently away from the vacuum interface.
 9. Put the ion source on a clean, safe surface.

Clean the Ion Source Surfaces



WARNING! Hot Surface Hazard. Before any maintenance procedures are started, let the temperature of the OptiFlow Turbo V ion source decrease for at least 60 minutes. Some surfaces of the ion source and vacuum interface become hot during operation.



WARNING! Electrical Shock Hazard. Before this procedure is started, remove the ion source from the mass spectrometer. Obey all electrical safe work practices.

Prerequisite Procedures

- [Remove the Ion Source](#).
- [Remove the Probe](#).

Clean the surfaces of the ion source after a spill or when they become dirty.

- Wipe the surfaces of the ion source with a soft, damp, cloth.

Remove the Probe



WARNING! Hot Surface Hazard. Before any maintenance procedures are started, let the temperature of the OptiFlow Turbo V ion source decrease for at least 60 minutes. Some surfaces of the ion source and vacuum interface become hot during operation.



WARNING! Electrical Shock Hazard. Before this procedure is started, remove the ion source from the mass spectrometer. Obey all electrical safe work practices.

CAUTION: Potential System Damage. Do not let the protruding electrode touch any part of the ion source housing, to avoid damaging the electrode.

The probe can be removed quickly and easily, without tools.

Note: If the probe is not correctly installed in the ion source, then the high-voltage power for the mass spectrometer and the source exhaust system are turned off.

Prerequisite Procedures

- Remove the column from the probe. Refer to the section: [Install the Micro Column and Heater](#) or [Install the Nano Column Cartridge and Heater](#).
- [Remove the Ion Source](#).

1. If a Micro probe is used, then remove the top fitting, with the integrated PEEK ferrule, and the bottom fitting from the probe. Refer to the section: [Install the Electrode in a Micro Probe](#).
2. Remove the electrode from the probe, and then put it on a safe, clean surface. Refer to the sections: [Install the Electrode in a Micro Probe](#) or [Install the Electrode in the Nano Probe](#) or [Connect the Echo[®] MS+ System to the OptiFlow Turbo V Ion Source](#).

CAUTION: Potential System Damage: Make sure that the electrode is removed from the probe before removing the probe from the ion source. Otherwise the electrode tip can be damaged.

3. Loosen the knurled ring on the probe, and then gently pull the probe straight up out of the source housing.
4. Put the probe on a clean, safe surface.

Tip! During installation of the probe in the ion source, align the dot on the probe with the corresponding dot on the ion source housing.

Tip! Clean the probe while it is removed from the ion source. Wipe the surfaces with a poly swab or lint-free wipe soaked in methanol.

Clean the Electrode

CAUTION: Potential System Damage. Do regular tests of the LC back pressure to make sure that the electrode is not blocked. Factors that can cause more frequent blockages include sample type, mobile phase type, usage time, and liquid that collected and dried in the electrode. We recommend that the LC back pressure be tested with a new, clean electrode to set a baseline. Then, do regular tests and compare the results with the baseline. If the back pressure increases very much, then clean or replace the electrode.

Note: To clean the electrode for the Echo[®] MS system, refer to the section, "OPI Port Wash", in the document: *Echo[®] MS+ System User Guide*.

Ion Source Maintenance

Required Materials
<ul style="list-style-type: none">• LC-MS–grade methanol or LC-MS–grade isopropanol, at ambient temperature

1. Remove the probe, with the electrode, from the ion source.
2. Connect the probe to the LC system.
3. Use the LC system to flush the probe with methanol or isopropanol at a minimum flow rate of 1 mL/min, until the back pressure is stable.

Storage and Handling



WARNING! Environmental Hazard. Do not discard system components in municipal waste. To discard components correctly, obey local regulations.

The environmental requirements for the storage and transport of the ion source:

- Ambient temperature between $-30\text{ }^{\circ}\text{C}$ and $+60\text{ }^{\circ}\text{C}$ ($-22\text{ }^{\circ}\text{F}$ and $140\text{ }^{\circ}\text{F}$)
- Atmospheric pressure between 75 kPa to 101 kPa
- Relative humidity not exceeding 99%, non-condensing

Connect the Echo[®] MS+ System to the OptiFlow Turbo V Ion Source

A

Note: The ZenoTOF 7600+ system does not support the Echo[®] MS+ system.


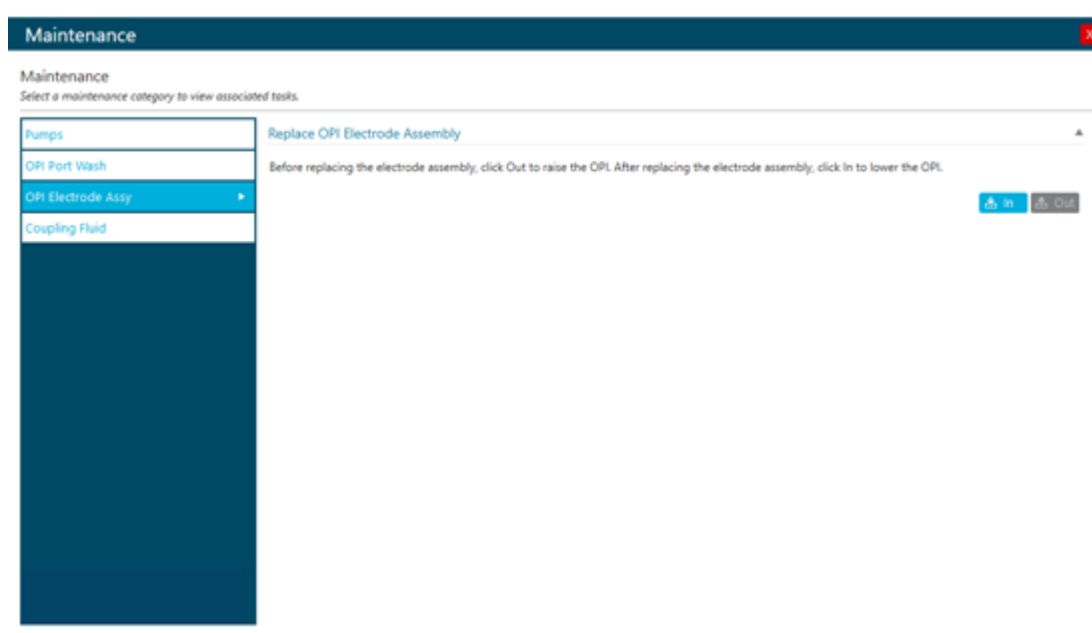
1. In the status panel, click  (**Direct device control**).
2. Click **Tools > Maintenance**.
3. In the Replace OPI Electrode Assembly section, click **Out**.

Figure A-1 Maintenance: OPI Electrode Assy



4. Deactivate the devices.

Connect the Echo® MS+ System to the OptiFlow Turbo V Ion Source

Figure A-2 Deactivate Devices



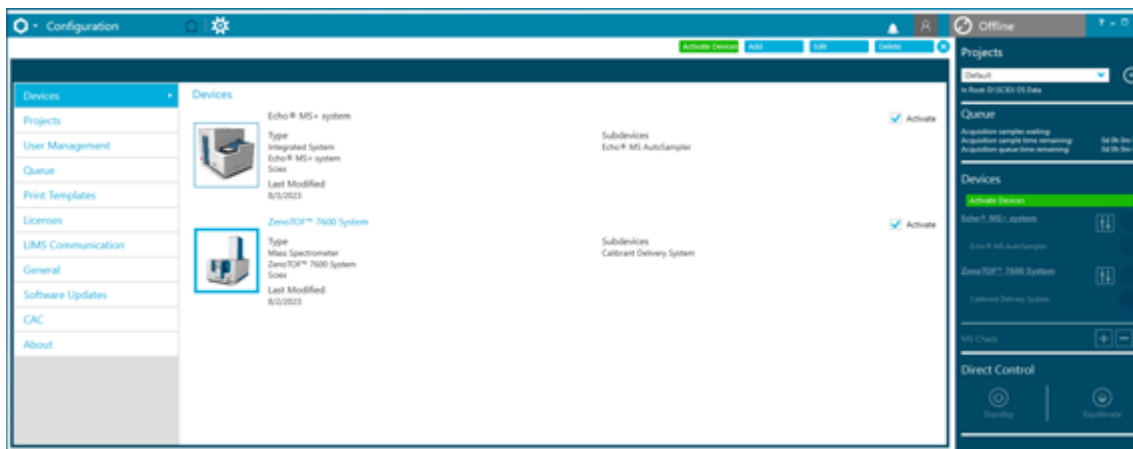


Figure A-3 Deactivate Devices

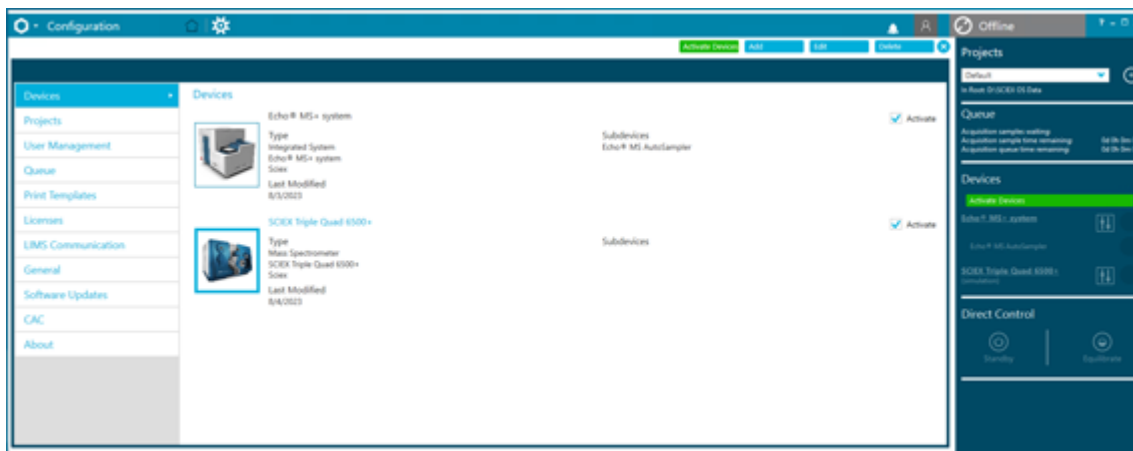

5. Disconnect the tubing from the Open Port Interface (OPI).


Figure A-4 Disconnect Tubing From the OPI

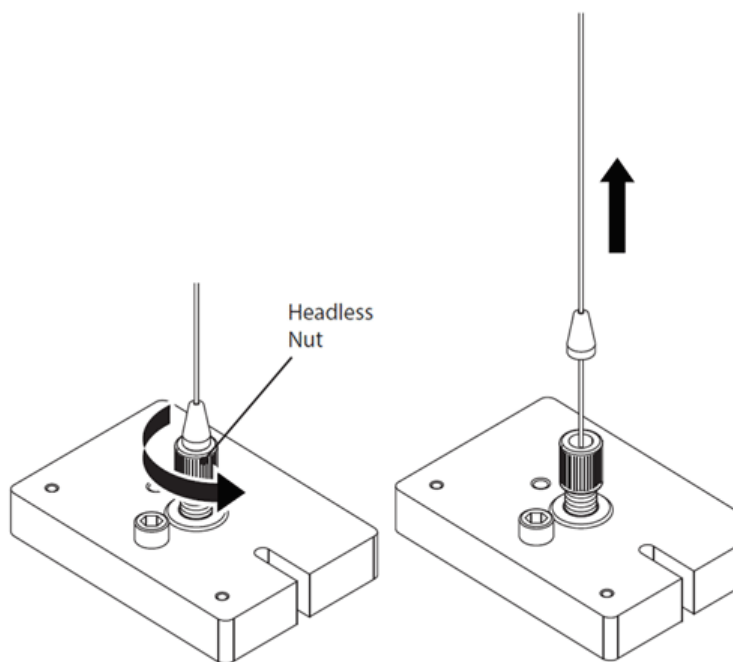

- a. Loosen the headless nut.
 - b. Disconnect the tubing from the OPI. Do not remove the headless nut.
6. Remove the probe from the OptiFlow Turbo V ion source.

Connect the Echo[®] MS+ System to the OptiFlow Turbo V Ion Source


Figure A-5 Remove the Probe

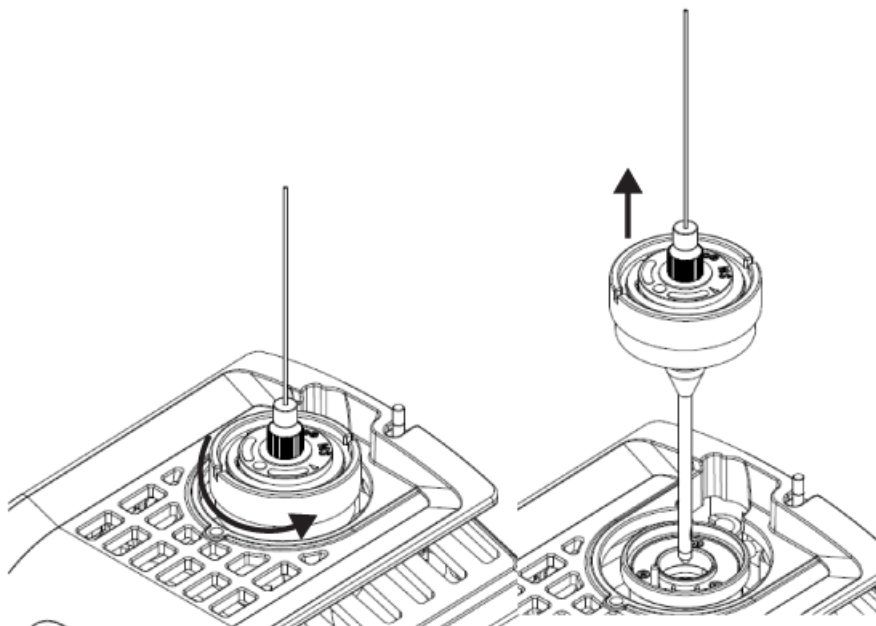

7. Remove the electrode assembly from the probe.

Note: Do not rotate the top fitting.

Figure A-6 Remove the Electrode Assembly

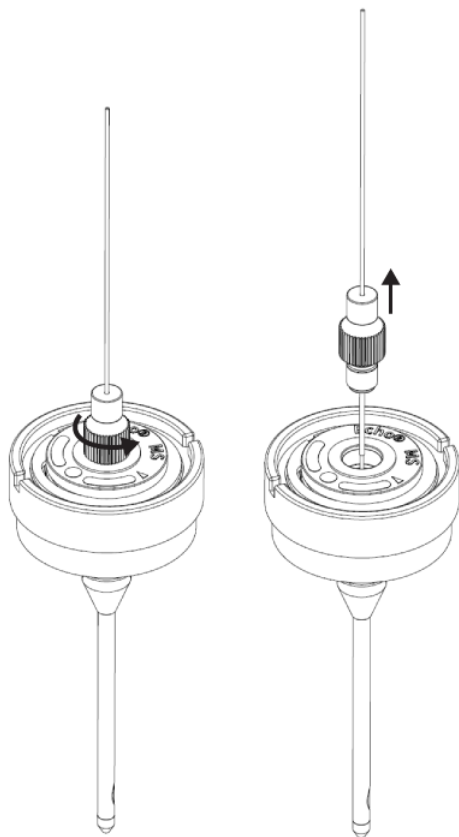


8. Remove the protective tubing from the new electrode assembly.

Note: Make sure to use the electrode assembly that comes with the probe.

9. Install and then fully tighten the new electrode assembly.

Connect the Echo[®] MS+ System to the OptiFlow Turbo V Ion Source


Figure A-7 Install the New Assembly

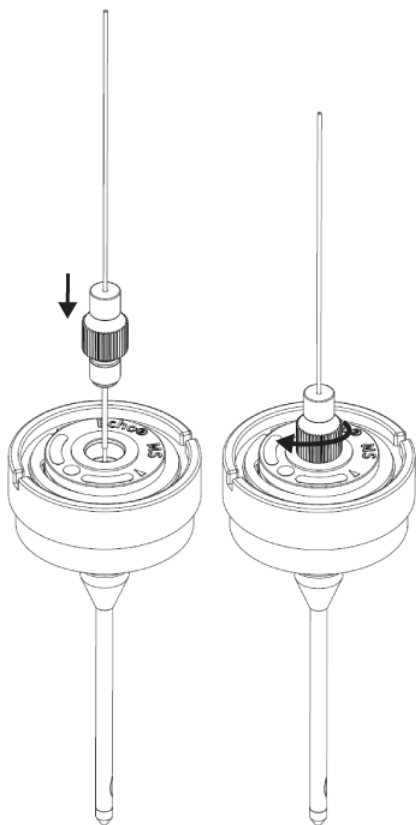

Note: The marking on the electrode assembly is not aligned with the marking on the probe.

Connect the Echo[®] MS+ System to the OptiFlow Turbo V Ion Source

Figure A-8 Electrode Assembly and Probe Markings Not Aligned

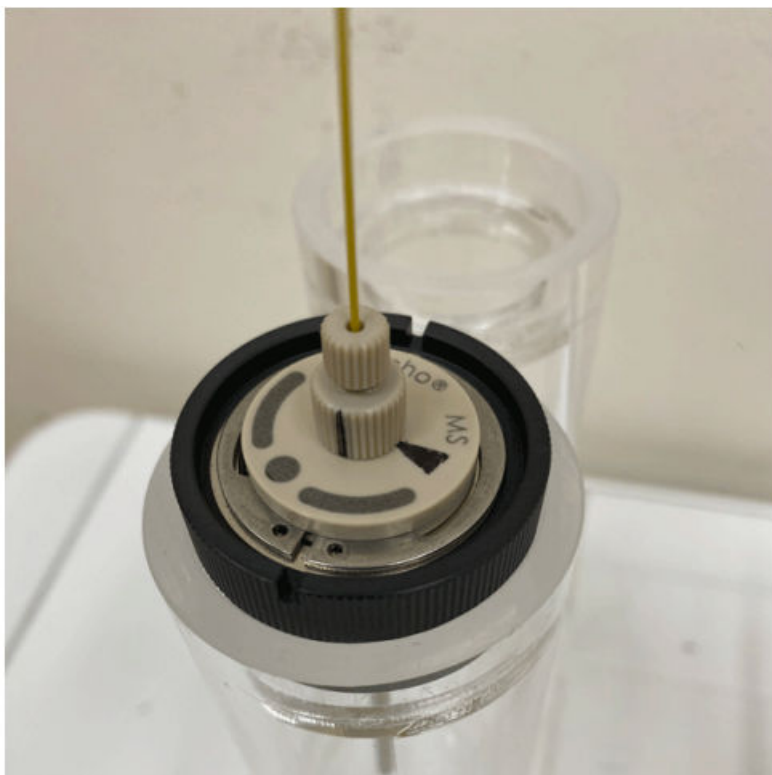

Note: Make sure that there is no gap between the bottom fitting and the probe.

Connect the Echo[®] MS+ System to the OptiFlow Turbo V Ion Source

10. Loosen the electrode assembly slowly, and then align the electrode assembly marking with the probe marking.


The marking identifies the position of the probe and electrode set after the assembly is optimized in the factory.

Figure A-9 Electrode Assembly and Probe Markings Aligned



11. Install the probe in the ion source.


Figure A-10 Install the Probe

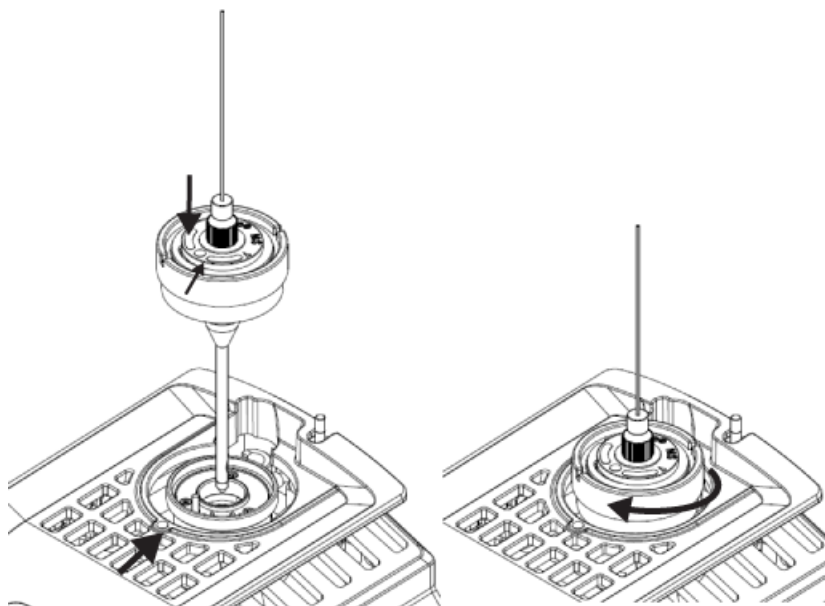

Connect the Echo[®] MS+ System to the OptiFlow Turbo V Ion Source

12. Connect the tubing to the OPI.

Figure A-11 Connect the Tubing

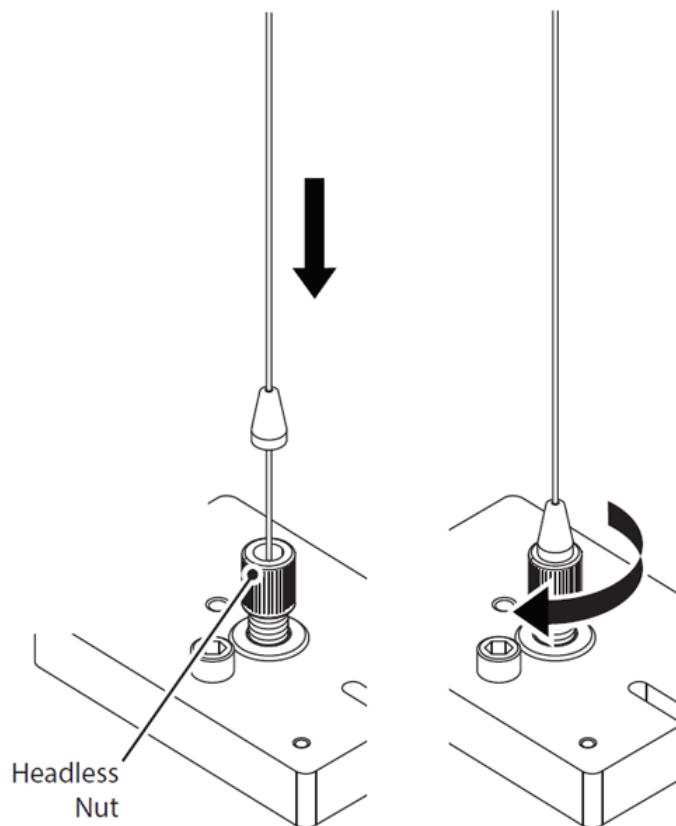

- Install the tubing through the headless nut until the stopper is flush with the nut.
- Tighten the nut.

Note: Make sure that there is no gap between the stopper and the headless nut.

13. Activate the devices.

Figure A-12 Activated Devices


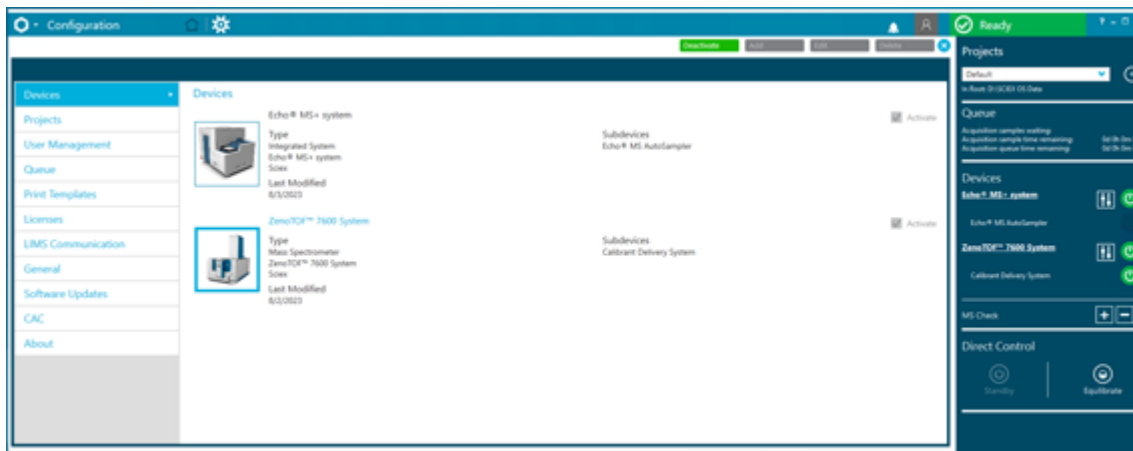


Figure A-13 Activated Devices

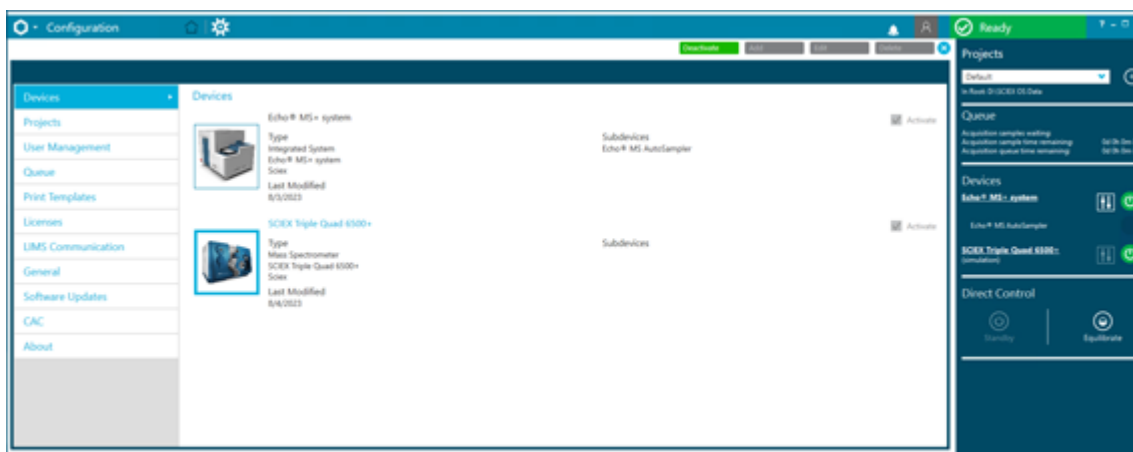


Note: An OPI Drip Sensor Fault detected error might occur if the system is initiated immediately after the electrode assembly is replaced. This error occurs when the carrier solvent flows out after the vacuum in the OPI is released. Wait at least 15 minutes for the sensor to dry and for the system to recover. To dry the OPI overflow sensor manually, refer to the section, "Clean the OPI Wash Tray" in the document: *Echo® MS+ System User Guide*.

14. In the Replace OPI Electrode Assembly section, click **In**.

Connect the Echo® MS+ System to the OptiFlow Turbo V Ion Source

Figure A-14 Maintenance: OPI Electrode Assy

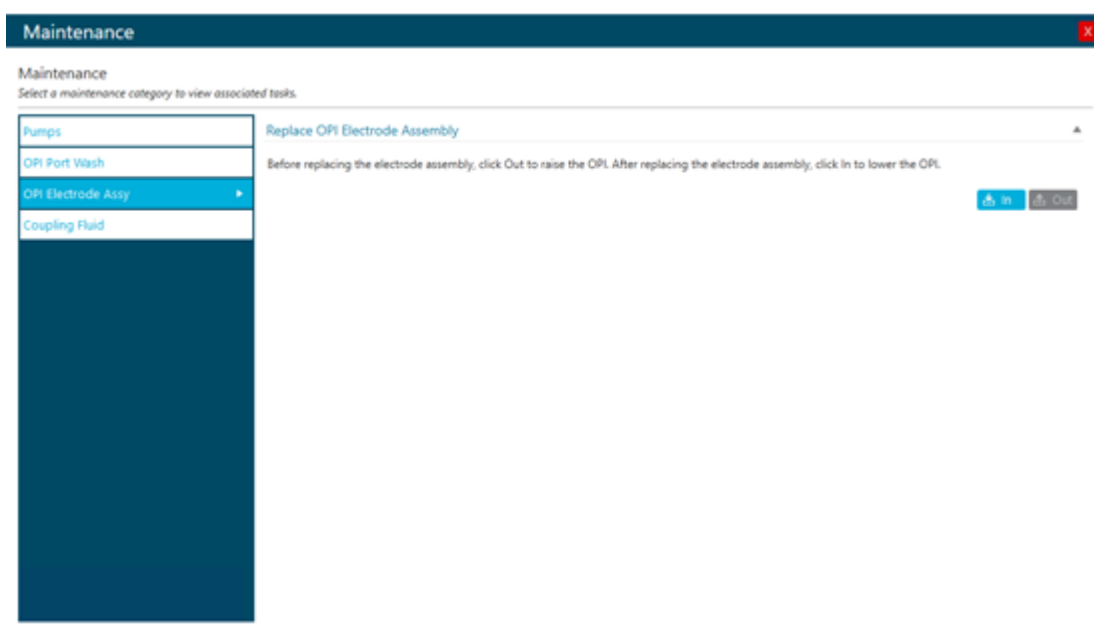




Labels on the Ion Source

B

In accordance with regulatory requirements, all warning labels shown on the ion source are documented in this guide. Warnings and labels on the ion source use international symbols.









Table B-1 Warning Labels

External Labels	Definition	Location
	ISO 7000-0434B (2004-1) CAUTION consult documentation	External
	CAUTION possibility of electric shock	External
	IEC 60417-5041 (2002-10) Caution hot surface	External









Glossary of Symbols

C











Note: Not all of the symbols in the following table are applicable to every instrument.

Symbol	Description
	Australian Regulatory Compliance Mark. Indicates that the product complies with Australian Communications Media Authority (ACMA) EMC and Electrical Safety Requirements.
	Alternating current
A	Amperes (current)
	Asphyxiation Hazard
	Authorized representative in the European community
	Biohazard
	CE Marking of Conformity
	cCSAus mark. Indicates electrical safety certification for Canada and USA.
	Catalog number





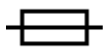





Glossary of Symbols

Symbol	Description
	Caution. Consult the instructions for information about a possible hazard. Note: In SCIEX documentation, this symbol identifies a personal injury hazard.
	China RoHS Caution Label. The electronic information product contains certain toxic or hazardous substances. The center number is the Environmentally Friendly Use Period (EFUP) date, and indicates the number of calendar years the product can be in operation. Upon the expiration of the EFUP, the product must be immediately recycled. The circling arrows show the product is recyclable. The date code on the label or product indicates the date of manufacture.
	China RoHS logo. The device does not contain toxic and hazardous substances or elements above the maximum concentration values and the device is an environmentally-friendly product that can be recycled and reused.
	Consult instructions for use.
	Crushing Hazard
	cTUVus mark for TUV Rheinland of North America
	Data Matrix symbol that can be scanned by a barcode reader to obtain a unique device identifier (UDI)
	Environmental Hazard






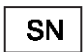

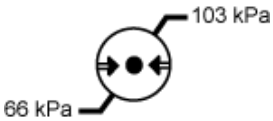
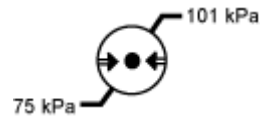
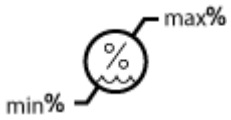

Glossary of Symbols

Symbol	Description
	Ethernet connection
	Explosion Hazard
	Eye Injury Hazard
	Fire Hazard
	Flammable Chemical Hazard
	Fragile
	Fuse
Hz	Hertz
	International safety symbol Caution, risk of electric shock (ISO 3864), also known as High Voltage symbol If the main cover must be removed, then contact a SCIEX representative to prevent electric shock.
	Hot Surface Hazard
	In Vitro Diagnostic Device








Glossary of Symbols

Symbol	Description
	Ionizing Radiation Hazard
	Keep dry. Do not expose to rain. Relative humidity must not exceed 99%.
	Keep upright.
	Lacerate/Sever Hazard
	Laser Radiation Hazard
	Lifting Hazard
	Magnetic Hazard
	Manufacturer
	Moving Parts Hazard
	Pacemaker Hazard. No access to people with pacemakers.

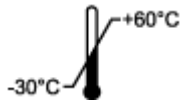Glossary of Symbols

Symbol	Description
	Pinching Hazard
	Pressurized Gas Hazard
	Protective Earth (ground)
	Puncture Hazard
	Reactive Chemical Hazard
	Serial number
	Toxic Chemical Hazard
	Transport and store the system within 66 kPa to 103 kPa.
	Transport and store the system within 75 kPa to 101 kPa.
	Transport and store the system within the specified minimum (min) and maximum (max) levels of relative humidity, noncondensing.
	Transport and store the system within $-30\text{ }^{\circ}\text{C}$ to $+45\text{ }^{\circ}\text{C}$.

Glossary of Symbols

Symbol	Description
	Transport and store the system within $-30\text{ }^{\circ}\text{C}$ to $+60\text{ }^{\circ}\text{C}$.
	USB 2.0 connection
	USB 3.0 connection
	Ultraviolet Radiation Hazard
	United Kingdom Conformity Assessment Mark
UKRP	United Kingdom Responsible Person
VA	Volt Ampere (apparent power)
V	Volts (voltage)
	WEEE. Do not dispose of equipment as unsorted municipal waste. Environmental Hazard
W	Watts (power)
	<i>yyyy-mm-dd</i> Date of manufacture

Contact Us

Customer Training

- In North America: NA.CustomerTraining@sciex.com
- In Europe: Europe.CustomerTraining@sciex.com
- Outside the EU and North America, visit sciex.com/education for contact information.

Online Learning Center

- [SCIEX Now Learning Hub](#)

SCIEX Support

SCIEX and its representatives maintain a staff of fully-trained service and technical specialists located throughout the world. They can answer questions about the system or any technical issues that might arise. For more information, visit the SCIEX website at sciex.com or contact us in one of the following ways:

- sciex.com/contact-us
- sciex.com/request-support

Cybersecurity

For the latest guidance on cybersecurity for SCIEX products, visit sciex.com/productsecurity.

Documentation

This version of the document supersedes all previous versions of this document.

To see this document electronically, Adobe Acrobat Reader is required. To download the latest version, go to <https://get.adobe.com/reader>.

To find software product documentation, refer to the release notes or software installation guide that comes with the software.

To find hardware product documentation, refer to the documentation that comes with the system or component.

The latest versions of the documentation are available on the SCIEX website, at sciex.com/customer-documents.

Note: To request a free, printed version of this document, contact sciex.com/contact-us.
